

INSTRUMENT INSTRUCTIONS

MICRO 500 CONDUCTIVITY METER MICRO 500 TDS METER



OPERATING INSTRUCTIONS

PT 144 MICRO 500 CONDUCTIVITY METER
(Conductivity/Temperature Meter)

PT 142 MICRO 500 TDS METER
(TDS /Temperature Meter)

TABLE OF CONTENTS

1	INTRODUCTION.....	4
2	DISPLAY AND KEYPAD FUNCTIONS.....	4
2.1	Display	4
2.2	Keypad.....	5
3	PREPARATION.....	6
3.1	Inserting and Removing Rubber Boot.....	6
3.2	Inserting the Batteries	7
3.3	Battery Replacement	7
3.4	Conductivity Electrode Information	8
3.5	Connecting the Probe to Meter	9
3.6	Switching the Meter On	10
3.7	Change Measurement Mode	11
4	CALIBRATION	12
4.1	Important Information on Meter Calibration	12
4.2	Preparing the Meter for Calibration.....	13
4.3	Selection of Automatic or Manual Calibration (PT 144 ONLY).....	14
4.4	Automatic Calibration (Conductivity).....	15
4.5	Manual Calibration (Conductivity or TDS).....	17
4.6	Temperature Calibration	18
5	MEASUREMENT.....	18
5.1	With Automatic Temperature Compensation (ATC)	18
5.2	Without ATC (Manual Temperature Compensation)	19
5.3	Taking Measurements	19
5.4	Using Manual Ranging Function.....	20
5.5	HOLD Function	21

6	ADVANCED SET-UP FUNCTIONS.....	22
6.1	Advanced Set-Up Overview.....	22
6.2	Select Cell Constant	24
6.3	Automatic Calibration (for PT 144).....	25
6.4	Setting the TDS Factory (for PT 142)	26
6.5	Temperature Coefficient	26
6.6	Normalization Temperature	27
6.7	Single-Point Calibration	27
6.8	Restore Factory Default Values.....	28
7	PROBE CARE AND MAINTENANCE.....	29
8	TROUBLE-SHOOTING GUIDE	30
9	ERROR MESSAGES	31
10	SPECIFICATIONS.....	32
11	ACCESSORIES.....	33

1 INTRODUCTION

Thank you for purchasing the Micro 500 TDS/Conductivity Meter. These microprocessor-based handheld meters deliver up to $\pm 0.5\%$ full-scale accuracy. It has a large custom LCD (liquid crystal display) for clear and easy reading.

PT 144 measures Conductivity ($\mu\text{S}/\text{mS}$) and Temperature ($^{\circ}\text{C}$) while the PT 142 measures Total Dissolved Solids (TDS) and Temperature ($^{\circ}\text{C}$). This sturdy meter measures up to 5 different ranges with auto-ranging capability that switches to appropriate measuring range automatically.

Your meter includes a conductivity electrode (cell constant $K = 1.0$) with built-in temperature sensor, a rubber boot, 4 alkaline 'AAA' batteries, instruction manual and soft pouch.


Please read this manual thoroughly before operating your meter.

To order other accessories and buffer standard solutions, please refer to the Accessories Section for more information.

2 DISPLAY AND KEYPAD FUNCTIONS

2.1 Display

The meter has a large custom LCD that consists of 4-digit segments and operation annunciators for $\mu\text{S}/\text{mS}$ (or ppm/ppt for PT 144 meter) and $^{\circ}\text{C}$ (Temperature). Other enunciators include 'HO' (when the HOLD function is activated) and 'LO' (low battery condition). See Figure below :-

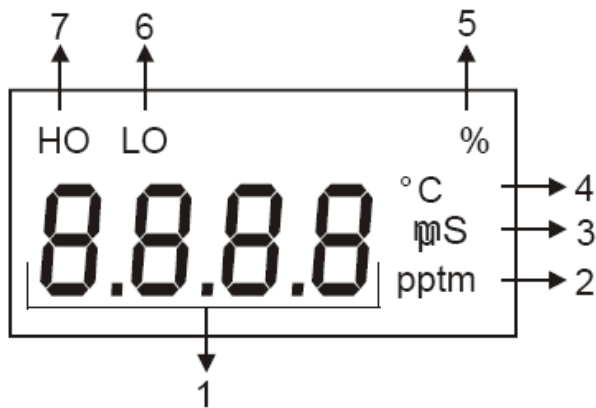

LCD and Customised Enunciators

- | | | | |
|---|-------------------------------------------------------------------------------------------------|---|--------------------------------------------------|
| 1 | Primary Display | 4 | Temperature indicator |
| 2 | Parts Per Million (ppm) or Parts Per Thousand (ppt) indicator – available for PT 144 meter only | 5 | Percentage indicator for Temperature Coefficient |
| 3 | Milli-Siemens/cm (mS) or micro-Siemens/cm (µS) indicator – available for PT 142 meter only. | 6 | Low Battery indicator |
| | | 7 | Hold (freeze) reading indicator. |

2.2 Keypad

The meter has six keys on its splash-proof keypad; ON/OFF, HOLD/ENTER, CAL, MODE, ▲ and ▼ keys. Some buttons have several functions depending on its mode of operation :-



Powers on and shuts off the meter. Takes you directly into measurement mode when meter is switched on.



Enters into calibration mode for Conductivity/TDS and Temperature.

To abort calibration or setup mode without confirming any set value.



HOLD: Freezes the measured reading. To activate, press HOLD key while in measurement mode. To release, press HOLD key again.

ENTER: Press to confirm values in calibration mode, and to confirm selections in SETUP mode.



In Calibration Mode: Press to scroll through calibration values.



In Setup Mode: Press to scroll through the setup sub-group programs.

Press ▲ key during conductivity measurement mode to activate manual ranging function. Each key press will move up higher conductivity range.

Selects measurement mode for conductivity/TDS and Temperature.



When pressed together with ON/OFF key, it will take you into the SETUP mode. This allows you to customize meter preferences such as selecting electrode's cell constant, normalization temperature, temperature coefficient factor, TDS factor (PT 144 only), automatic (PT 142 only) or manual calibration, single-point or multi-point calibrations, and to reset meter to factory default.

3 PREPARATION

3.1 Inserting and Removing Rubber Boot

(Rubber Boot sold separately)

- 1) To remove meter from rubber boot, push out from the bottom edges of meter until it is completely out of the boot. Ensure that the cables of conductivity electrode or temperature probe are not connected.




- 2) To insert meter into rubber boot, slide in from the top of meter before pushing the bottom edges of meter down to set it into position. Lift up the stand at the back of meter for bench top applications if necessary.

3.2 Inserting the Batteries


The battery compartment is found at the back of instrument as shown. To open the battery compartment :-

- 1) Push in the direction of arrow and lift up the cover.
- 2) Note the polarity of battery before inserting into position.
- 3) After replacement, place cover back and press down until it locks tight.

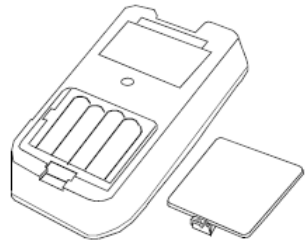

3.3 Battery Replacement

A 'LO' enunciator in the LCD alerts you when battery power is running low. See Figure below. Replace with the same type as recommended by the manufacturer :-

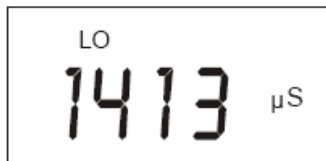

LO Battery Condition

Caution: Power off the meter when changing battery.

3.4 Conductivity Electrode Information

The meter is supplied with a conductivity/TDS electrode with a BNC connector. This conductivity/TDS electrode comes with Stainless Steel rings, cell constant of $K = 1.0$, and a built-in temperature sensor for Automatic Temperature Compensation (ATC). Its specially designed Ultem-body housing has good chemical resistant properties. It provides fast temperature response and reduces air entrapment, which makes it easy to obtain accurate, stable readings.

The probe materials used which have good chemical durability include :

- 1 Polyetherimide (Ultem) – protective probe guard
- 2 Polybutylterphalate (Valox) – sensor housing
- 3 Stainless Steel (SS 304) – 2 steel bands


Proper use of probe is essential to ensure that the optimum measurement is taken in a short time.

The removable protective plastic probe guard is meant for simple periodic maintenance and it must be kept intact during measurement and calibration.


Always immerse the probe beyond upper steel band.

NOTE :

- 1 DO NOT remove the protective probe guard during measurement and calibration as it may affect your readings.

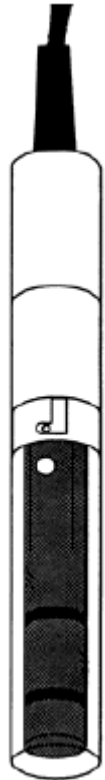

- 2 We recommend that you do not submerge the probe above the protective guard. You can submerge the cable for brief periods of time, but not continuously :-

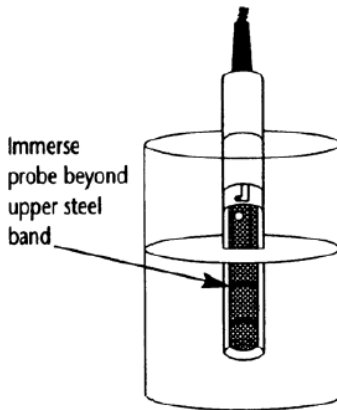


See *Section 7 – ‘Probe Care and Maintenance’* for more information.

3.5 Connecting the Probe to Meter

- 1) To connect electrode into meter, align the connector slots with the posts of meter’s socket and rotate connector clockwise until it locks.
- 2) To remove, simply rotate the connector in anti-clockwise direction until it unlocks, and slide the connector off the socket.
- 3) Insert the mini phono-jack of temperature sensor into the socket on the meter as shown below.
- 4) Unplug the phono-jack when not in use or when you want to measure Conductivity or TDS without any temperature compensation (Manual Temperature Compensation, see Section 5.2).

CAUTION: Do not pull or force on the probe cord or the probe wires might disconnect.

NOTE: Keep connectors clean.
Do not touch connector with soiled hands.

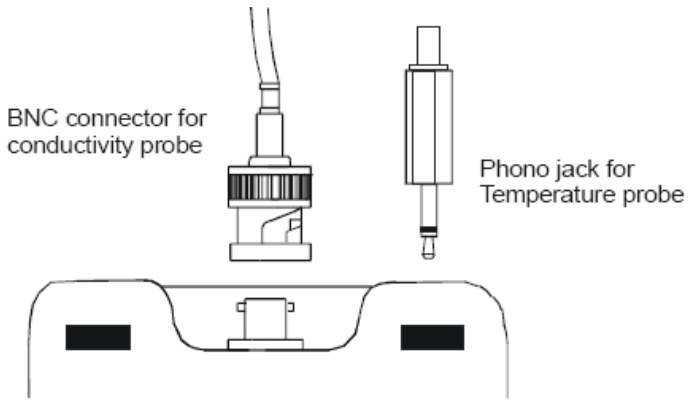


Connection for Conductivity and Temperature Probes


3.6 Switching the Meter On

When switching the meter on, it will go through a series of display, showing the various setup parameters.

For Micro 500 Conductivity Meter :-

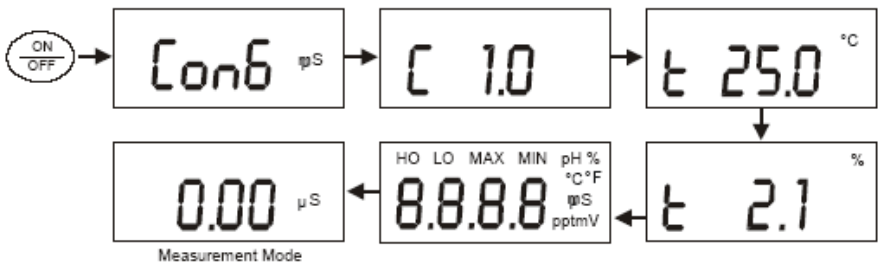

For Micro 500 TDS Meter :-

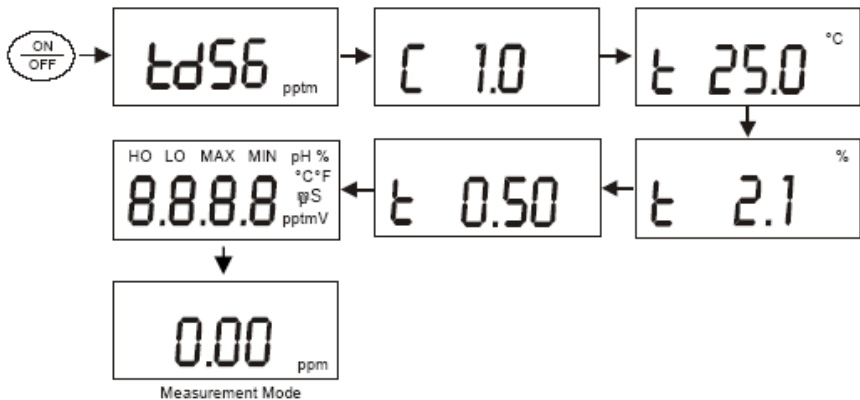


Press **ON/OFF** key to power up your meter.

- 1) The first screen shows [**Con 6**] (or [**tdS 6**]) depending on the meter type.
- 2) Second screen shows [**C 1.0**] which is the conductivity cell constant, k. You can select different cell constant of 0.1, 1.0 or 10.0. Refer to Section on Advanced Set-Up. Default value is $k = 1.0$.
- 3) Third screen shows [**t 25.0 °C**] which is the Normalization Temperature. You can set Normalization Temperature at either 25°C or 20°C. Refer to Section on Advance Setup. Default value is 25°C.
- 4) Fourth screen shows [**t 2.1%**] which is the Temperature Coefficient. You can customise the meter with different Temperature Coefficient value from 0.0 to 3.0 %/°C from the Advance Setup mode. Default value is 2.1% / °C.
- 5) All LCD segments will light up for two seconds, and change into measurement mode.
- 6) You are now ready for conductivity measurement.

3.7 Change Measurement Mode

To switch between Conductivity/TDS measurement mode and Temperature measurement mode, simply press the **MODE** key.

The customized enunciator helps indicate the measurement parameter you are in.

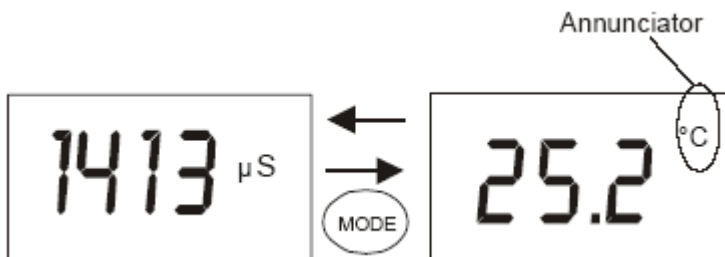

4 CALIBRATION

4.1 Important Information on Meter Calibration

Your meter has five measuring ranges. You can calibrate one point in each of the measuring ranges (up to five points). If you are measuring values in more than one range, make sure to calibrate each of the ranges you are measuring.

The following table lists the corresponding conductivity and TDS ranges. You should calibrate each range using a solution that falls between the values in the 'recommended calibration solution range' column :-

Conductivity Range	Recommended Calibration Solution Range	TDS Range	Recommended Calibration Solution Range
0.00 \rightarrow 20.00 μS	6.00 to 17.00 μS	0.00 \rightarrow 10.00 ppm	3.00 to 8.50 ppm
0.0 \rightarrow 200.0 μS	60.0 to 170.0 μS	10.0 \rightarrow 100.0 ppm	30.0 to 85.0 ppm
0 \rightarrow 2000 μS	600 to 1700 μS	100 \rightarrow 1000 ppm	300 to 850 ppm
0.00 \rightarrow 20.00 mS	6.00 to 17.00 mS	1.0 \rightarrow 10.00 ppt	3.00 to 8.50 ppt
0.0 \rightarrow 200.0 mS	60.0 to 170.0 mS	10.0 \rightarrow 200 ppt	30.0 to 170 ppt

When you recalibrate your meter, old calibrations are replaced on a range basis. For example, if you previously calibrated your conductivity meter at 1413 μS in the 0 to 2000 μS range and you recalibrate at 1500 μS (also in the 0 to 2000 μS range), the meter will replace the old calibration data (1413 μS) in that range. The meter will retain all calibration data in other ranges.

To completely recalibrate your meter, or when you use a replacement probe, it is best to clear all calibration data. To erase all the old conductivity or TDS calibration data completely, see *Section 6.8 – Restore Factory Default Values*.

4.2 Preparing the Meter for Calibration

Before starting calibration, make sure you are in the correct measurement mode.

For best results, select a standard value close to the sample value you are measuring. Alternatively use a calibration solution value that is approximately 2/3 the full-scale value of the measurement range you plan to use. For example, in the 0 to 2000 μS conductivity range, use a 1413 μS solution for calibration.

Calibrate to all measurement ranges to ensure the highest accuracy throughout all measurement range. Note that the meter will not accept calibration values less than 40 $\mu\text{S}/\text{cm}$ (20 ppm). All new calibration values will automatically override existing data.

If you are measuring in solutions with Conductivity lower than 100 $\mu\text{S}/\text{cm}$ or TDS lower than 50 ppm, calibrate the meter at least once a week to get good accuracy. If you are measuring in the mid ranges and you wash the probe in deionized water and store it dry, calibrate the meter once a month. If you take measurements at extreme temperatures, calibrate at least once a week.

Ensure that you use new Conductivity standard solutions or sachets during calibration. Do not reuse standard solutions as it may be contaminated and affect the calibration and accuracy of measurements. Use fresh calibration solution each time you calibrate your meter. Store solutions in a dry and cool environment if possible.

Always rinse the probe with either deionized water or rinse solution before and after each calibration/sample measurement to avoid cross-contamination. For details please refer to *Section 7 - Probe Care and Maintenance*.

NOTE: These meters are factory set to a temperature coefficient of 2.1% / °C. For most applications this will provide good results. To set the temperature coefficient to different value, see Section 6.5 – Temperature Coefficient. Also, see Addendum 3 - Calculating Temperature Coefficient to determine the appropriate temperature coefficient for your solution.

NOTE: The factory default value for normalization temperature is 25 °C. If you need to normalize to a value other than 25 °C, see Section 6.6 – Normalization Temperature.

4.3 Selection of Automatic or Manual Calibration (PT 144 ONLY)

The Micro 500 Conductivity meter is capable of performing either automatic or manual calibration.

In the automatic calibration mode, the meter (PT 144 only) automatically detects and verifies the appropriate known calibration standards solutions being calibrated before accepting these particular calibration standards as one of its calibration values in a specific measurement range. This automatic calibration mode frees you from cumbersome calibration procedure.

The known calibration standards used for automatic calibration are :-

Meter	Normalisation Temperature	Calibration Standards (Range)
PT 142	25 °C	1) 84 µS (for 0 – 200 µS/cm) 2) 1413 µS (for 0 – 2000 µS/cm) 3) 12.88 mS (for 0.00 – 20.00 mS/cm) 4) 111.8 mS (for 0.0 – 200.0 mS/cm)
	20 °C	1) 76 µS (for 0 – 200 µS/cm) 2) 1278 µS (for 0 – 2000 µS/cm) 3) 11.67 mS (for 0.00 – 20.00 mS/cm) 4) 102.1 mS (for 0.0 – 200.0 mS/cm)

Table 1: Conductivity Calibration Standards for Auto Calibration

In the manual calibration, non-standard calibration values can be used for calibration. You can manually input the appropriate values as your desired calibration standards in each specific range. This is useful when you have a customized calibration standard specifically unique for your application.

To select Automatic or Manual Calibration settings, refer to *Section 6.3 – Automatic Calibration* for more information.

4.4 Automatic Calibration (Conductivity)

In the Automatic Calibration mode, the meter is capable of accepting either single-point or up to four points for multi-point calibration with maximum of one point per specific measurement range. For the known calibration standard values refer to *Table 1 in Section 4.3*.


- 1) If necessary, press **MODE** key to select conductivity mode.
- 2) Rinse the probe thoroughly with deionized water or a rinse solution, then rinse with a small amount of calibration standard.

NOTE: For Automatic Calibration you must use one of the calibration standards listed in Table 1.

- 3) Dip the probe into the calibration standard. Immerse the probe tip beyond the upper steel band (see *Figure in Section 3.4*). Stir the probe gently to create a homogeneous sample. Allow time for the reading to stabilize.
- 4) Press **CAL** key to enter conductivity calibration mode. The **[CA]** indicator will appear for 1½ seconds, and a value will appear flashing.

NOTE: To exit calibration without confirmation, press CAL key again to go back to measurement mode.

- 5) Wait for the value to stabilize and press **ENTER** key. The calibration standard value will appear for three seconds. If the calibration is successfully performed, a **[donE]** will be displayed for about three seconds, and the meter returns to measurement mode.

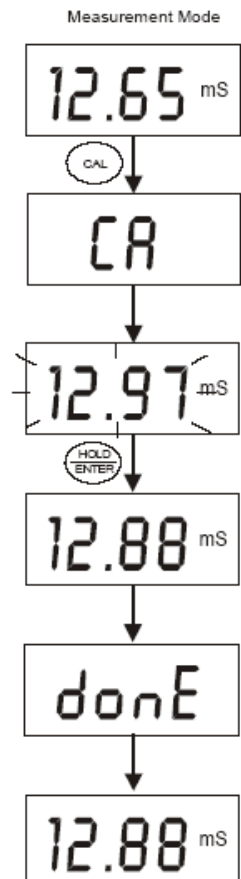

- 6) To perform the next point calibration in the multi-point calibration, repeat steps 1 - 5 again until all points have been calibrated if necessary.

IMPORTANT NOTES :

- 1 *Meter allows a tolerance range of $\pm 40\%$ of its calibration standard. An error message 'Err 1' will be displayed for three seconds if you attempt to calibrate with a solution whose value is outside the tolerance range. For example: for 1413 μS conductivity calibration standard, 40% tolerance is from 848 μS to 1978 μS .*
- 2 *If the temperature ($t^{\circ}\text{C}$) of the conductivity calibration solution is below 0°C or above 50°C ($0^{\circ}\text{C} < t^{\circ}\text{C} > 50^{\circ}\text{C}$), an error message 'Err 2' will be displayed when performing the auto calibration, and meter will return to measurement mode.*
- 3 *All new calibration data will over-ride existing stored calibration data for each measuring range calibrated.*
- 4 *It is important to use new conductivity calibration standards.*
- 5 *Low conductivity standard solution (less than 20 $\mu\text{S}/\text{cm}$) cannot be available easily. Such low conductivity standard will be contaminated as soon as it is exposed to the air therefore exercise caution during calibration in the first measurement range (0.00 to 20.0 $\mu\text{S}/\text{cm}$).*

4.5 Manual Calibration (Conductivity or TDS)

In Manual Calibration mode, you can use customized conductivity calibration standards (specific to your own application) and calibrate the meter. The following example shows calibration sequence to 12.00 mS conductivity calibration standard.

The procedure is similar for PT 144 and PT 142 meters :-

- 1) If necessary, press **MODE** key to select conductivity mode.
- 2) Rinse the probe thoroughly with de-ionized water or a rinse solution, then rinse with a small amount of calibration standard.
- 3) Dip the probe into the calibration standard. Immerse the probe tip beyond the upper steel band (see Figure in Section 3.4). Stir the probe gently to create a homogeneous sample. Allow time for the reading to stabilize.
- 4) Press **CAL** key to enter conductivity calibration mode. The [CA] indicator will appear for 1½ seconds, and a value will appear flashing.

NOTE: To exit calibration without confirmation, press CAL key again to go back to measurement mode.

- 5) Wait for the value to stabilize and press ▲ or ▼ key and adjust the value to the calibration standard used.
- 6) Press the **ENTER** key. The [CO] indicator will appear for 1½ seconds, and the calibration is successfully performed. The meter returns to measurement mode.
- 7) To perform the next point calibration in the multi-point calibration for next range, repeat steps 1 - 6 again until all points have been calibrated if necessary.



4.6 Temperature Calibration

The Conductivity electrode (Refer to Accessories section for order number) has a built-in temperature sensor for ATC. The temperature sensor is factory calibrated to the meter. Calibrate your sensor only if you suspect temperature errors that may have occurred over a long period of time or if you have a replacement probe.

- 1) Make sure that the phono-jack (for temperature measurement) is properly connected to the meter. *See Figure in Section 3.5.*
- 2) Switch on the meter and if necessary, press **MODE** key to select temperature measurement mode. *See Section 3.7.*
- 3) Press **CAL** key to start temperature calibration process.
- 4) Dip the probe into a solution with known temperature (for example, a temperature bath). Allow time for the temperature to stabilise.
- 5) Wait for the value to stabilize and press **▲** or **▼** key and adjust the value to the solution temperature.
- 6) Press the ENTER key. The [CO] indicator will appear for 1½ seconds, and the reading will stop flashing. The temperature calibration is successfully performed. The meter returns to measurement mode.

NOTE: To exit calibration without confirmation, press CAL key again to go back to measurement mode.

NOTE: You can offset the temperature reading up to $\pm 5^{\circ}\text{C}$ from the original (default) reading.

5 MEASUREMENT

The meter is capable of taking measurements with automatic or manual temperature compensation.

5.1 With Automatic Temperature Compensation (ATC)

For ATC, make sure the phono-jack of the probe (*see Figure in Section 3.5*) is securely inserted.


The conductivity/TDS reading displayed will be compensated for according to the normalization temperature (20°C or 25°C) selected. *See Section 6.6 – Normalization Temperature.*

5.2 Without ATC (Manual Temperature Compensation)

For manual temperature compensation, simply unplug the probe's phono-jack (not BNC) from the meter.

To use manual temperature compensation, you need to enter the temperature value of your process into the meter. This is the value at which the reading will manually temperature compensate. You can select any temperature between 0 and 50°C (32 to 122°F). Default value is 25°C :-

- 1) Make sure that the phono-jack (for temperature measurement) is disconnected from the meter. *See Figure in Section 3.5.*
- 2) Switch on the meter and if necessary, press **MODE** key to select temperature measurement mode. *See Section 3.7.*
- 3) Press **CAL** key to start temperature calibration process.
- 4) The 'CA' will appear momentarily and a temperature value will start flashing.
- 5) Check the temperature of your sample using an accurate thermometer. Wait for the value to stabilize and press ▲ or ▼ key and adjust the value to the reference thermometer used.
- 6) Press the **ENTER** key. The [CO] indicator will appear for 1½ seconds, and the reading will stop flashing. The temperature calibration is successfully performed. The meter returns to measurement mode.

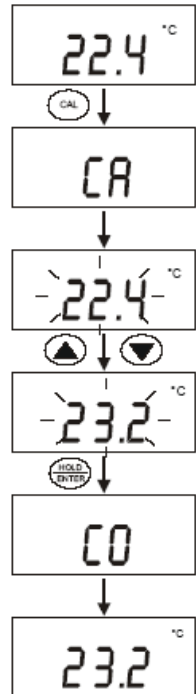

5.3 Taking Measurements

To take readings :-

- 1) Rinse the probe with deionized or distilled water before use to remove any impurities adhering to the probe body. Shake or air dry. To avoid contamination or dilution of your sample, rinse probe with a small volume of your sample liquid.
- 2) Press **ON** to switch on meter.

- 3) Dip the probe into the sample.
- 4) Allow time for the reading to stabilize. Note the reading on the display.


NOTE: When dipping the probe into the sample take care to ensure that the liquid level is above its upper steel band. Stir the probe gently in the sample to create a homogenous sample. See Figure in Section 3.4.

5.4 Using Manual Ranging Function

By default your meter has auto-ranging ability and would automatically selects the range in which your readings appear.

However, you may also manually select a specific range you want to measure. This is possible by simply pressing ▲ key successively for each measurement range. The five ranges are :-

Conductivity Range (PT 142)	TDS Range (PT 144) (if TDS factor is 0.5)
0 – 20.00 $\mu\text{S}/\text{cm}$	0 – 10.00 ppm
0 – 200.0 $\mu\text{S}/\text{cm}$	0 – 100.0 ppm
0 – 2000 $\mu\text{S}/\text{cm}$	0 – 1000 ppm
0 – 20.00 mS/cm	0 – 10.00 ppt
0 – 200.0 mS/cm	0 – 100 ppt

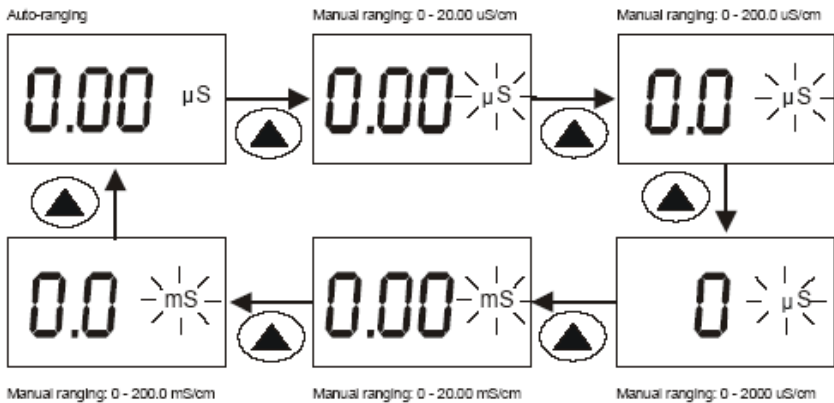

NOTE: If the value of the solution you are measuring is higher than the range selected [Or] will appear on the primary display. Press RANGE until the correct range is selected.

The meter resets to the Auto-ranging function once it is turned off. You will have to reset the manual ranging function each time you turn the meter off.

5.5 HOLD Function

This feature lets you freeze the display for a delayed observation. **HOLD** can be used any time in measurement mode.

- 1) To hold a measurement, press the **HOLD** key while in measurement mode. [HO] will appear on the display.
- 2) To release the held value, press the **HOLD** again. Continue to take measurements.



NOTE: This meter shuts off automatically after 20 minutes of non-use.





If the meter is shut off either automatically or manually, the HOLD value will be lost.

6 ADVANCED SET-UP FUNCTIONS

6.1 Advanced Set-Up Overview

The advanced setup mode lets you customize your meter's preferences and defaults. To enter advanced set-up mode :-

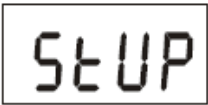- 1) Make sure that the meter is switched-off.
- 2) Press **ON** and **MODE** key simultaneously, holding both keys for two seconds. First release **ON** key first before releasing the **MODE** key.
- 3) [**StUP**] indicator will appear momentarily and [**CELC**] will appear next.
- 4) Overviews of Set-Up Menu are as follows :-

	Enter Set-Up Page
	Select Cell Constant Choice of k = 0.1, 1.0 and 10.0 Default value = 1.0
	PT 144 ONLY - Select 'Automatic Calibration'. 'Yes' for auto-calibration, and 'No' for manual calibration. Default value is 'Yes'.
	Adjust Temperature Coefficient value from 0.0 to 3.0%/°C Default value = 2.1%/°C

t.nr °C

Select 'Normalisation Temperature'
Choice of either 20°C or 25°C
Default value is 25°C

t.d.S

PT 142 ONLY - Adjust TDS factor from 0.4 to 1.0
Default value is 0.5


S.P.CA

Select 'Single Point Calibration'
Choice of 'Yes' or 'No'
Default value is 'Yes'


UrSt

User reset to factory defaults
Choice of 'Yes' or 'No'
Default value is 'No'

Overview of Advanced Set-Up

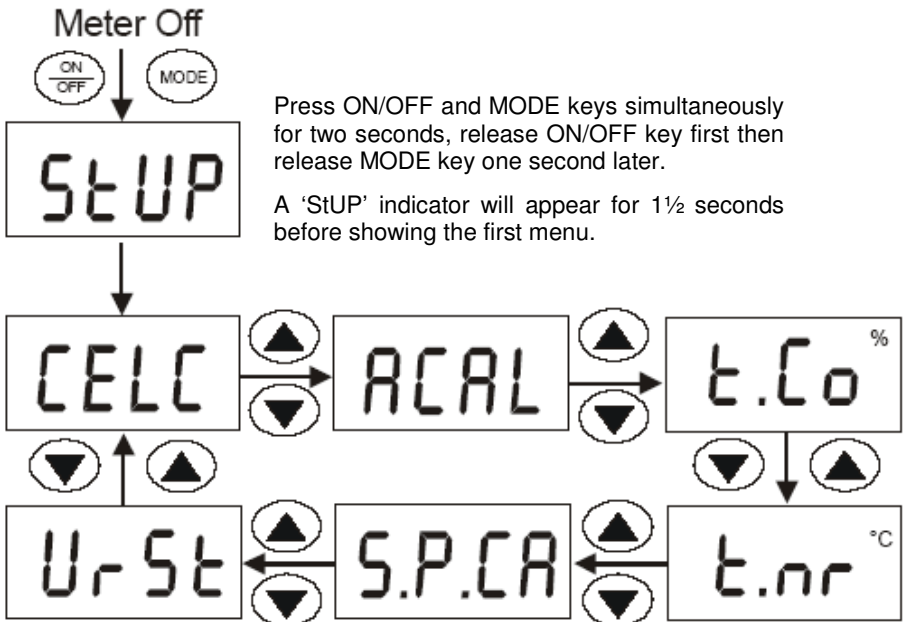


Overview of PT 144 Set-Up Menu

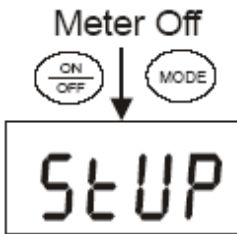

Press ON/OFF and MODE keys simultaneously for two seconds, release ON/OFF key first then release MODE key one second later.

A 'StUP' indicator will appear for 1½ seconds before showing the first menu.

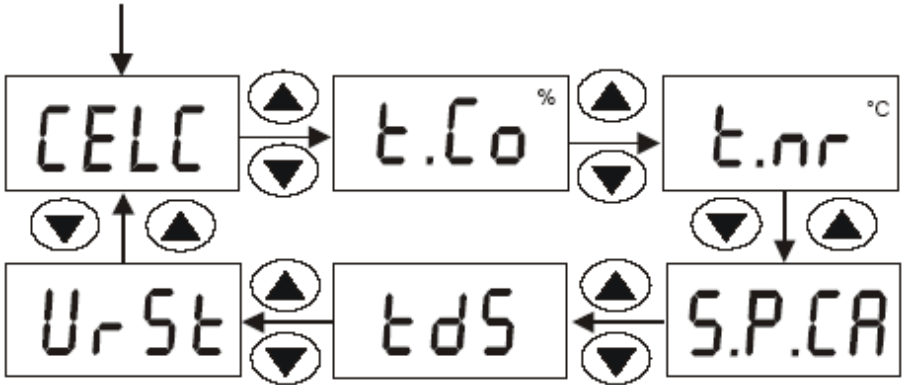

Overview of PT 142 Set-Up Menu


6.2 Select Cell Constant

This meter lets you select a cell constant of $K = 1.0, 10$ or 0.1 :-

- Use a cell of $K = 1.0$ for midrange measurements
- Use a cell of $K = 10$ for high range measurements (above 20 mS or 10 ppt).
- Use a cell of $K = 0.1$ for low range measurements (below 20 μ S or 10 ppm).

The cell included with your meter has a cell constant of $K = 1.0$.

- 1) Enter the advanced set-up as described in Section 6.1.
- 2) Press \blacktriangle or \blacktriangledown key until [CELC] appears on the LCD. Press **ENTER** key.

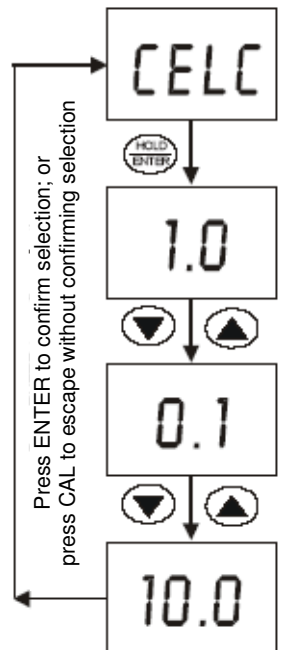


- 3) Press ▲ or ▼ key to select either '1.0', '0.1' or '10.0'. Ensure the cell constant selected correspond with the conductivity electrode you are using.
- 4) Press **ENTER** key to select. The meter will take you back to the menu, [**CELC**].
- 5) Press ▲ or ▼ key to move to the next menu or press **CAL** to exit to measurement mode.

6.3 Automatic Calibration (for PT 144)

The automatic calibration allows you to quickly calibrate the meter to any of the four widely used conductivity calibration standards. For a list of calibration standards refer to *Table 1 in Section 4.3*.

In the manual calibration mode, you can use you own customized conductivity calibration standard to calibrate this meter :-

- 1) Enter the advanced setup as described in Section 6.1.
- 2) Press ▲ or ▼ key until [**ACAL**] appears on the LCD. Press **ENTER** key.
- 3) Press ▲ or ▼ key to select either [**Yes**] or [**No**].
- 4) Press [ENTER] key to select. The meter will take you back to the menu, [**ACAL**].
- 5) Press ▲ or ▼ key to move to the next menu or press **CAL** to exit to measurement mode.

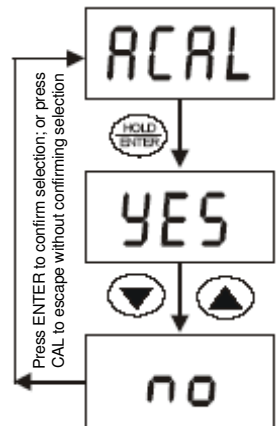

6.4 Setting the TDS Factory (for PT 142)


The concentration of salts dissolved in solution increases the conductivity of that solution. This relationship varies from salt to salt and is roughly linear over a given range for a given salt. The TDS conversion factor is the number used by the meter to convert from conductivity to TDS.

To calculate the TDS conversion factor refers to *Addendum 2 – Calculating TDS Conversion Factor*.

You can also look up at various Chemical reference books for TDS factor for various types of salt.

You can set the TDS conversion factor between 0.4 and 1.0; meter default is 0.5.

- 1) Enter the advanced setup as described in Section 6.1.
- 2) Press ▲ or ▼ key until [**tdS**] appears on the LCD. Press **ENTER** key.
- 3) Press ▲ or ▼ key to select a value between 0.4 to 1.0.
- 4) Press **ENTER** key to select. The meter will take you back to the menu, [**tdS**].
- 5) Press ▲ or ▼ key to move to the next menu or press **CAL** to exit to measurement mode.

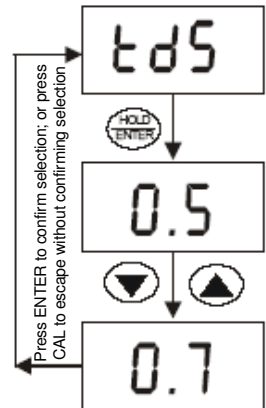


6.5 Temperature Coefficient

The temperature coefficient is the amount of change in conductivity per degree of temperature; it is expressed in percent per °C. Entering the exact temperature coefficient of your solution lets you accurately compensate temperature for almost any solution. You can adjust 0.0 to 3.0% per °C.

Meter default is 2.1% per °C.

- 6) Enter the advanced setup as described in Section 6.1.
- 7) Press ▲ or ▼ key until [t.Co %] appears on the LCD. Press **ENTER** key.
- 8) Press ▲ or ▼ key to select a value between 0.0 to 3.0.
- 9) Press **ENTER** key to select. The meter will take you back to the menu, [t.Co %].


Press ▲ or ▼ key to move to the next menu or press **CAL** to exit to measurement mode.

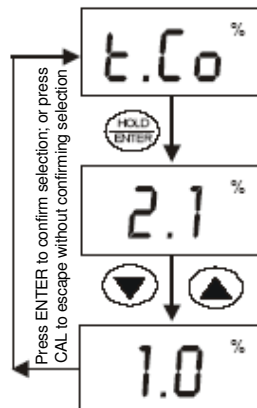

6.6 Normalization Temperature

You can set the meter to normalize its conductivity measurements to a standard temperature of either 25°C or 20°C. The default value is 25°C.

- 1) Enter the advanced setup as described in Section 6.1.
- 2) Press ▲ or ▼ key until [t.nr °C] appears on the LCD. Press **ENTER** key.
- 3) Press ▲ or ▼ key to select either [25.0 °C] or [20.0 °C].
- 4) Press **ENTER** key to select. The meter will take you back to the menu, [t.nr °C].
- 5) Press ▲ or ▼ key to move to the next menu or press **CAL** to exit to measurement mode.

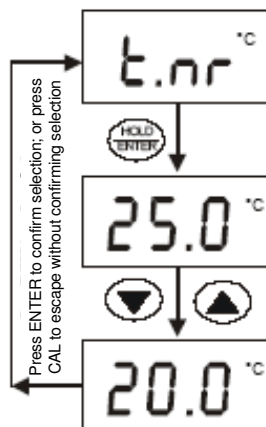


6.7 Single-Point Calibration

Single-point calibration refers to calibrating one conductivity value and uses it for the entire five conductivity ranges.

By selecting [no] to single-point calibration, you can perform calibration for each conductivity range :-

- 1) Enter the advanced setup as described in Section 6.1.
- 2) Press ▲ or ▼ key until [S.P.CA] appears on the LCD. Press **ENTER** key.
- 3) Press ▲ or ▼ key to select either [**Yes**] or [**no**].
- 4) Press **ENTER** key to select. The meter will take you back to the menu, [S.P.CA].
- 5) Press ▲ or ▼ key to move to the next menu or press **CAL** to exit to measurement mode.

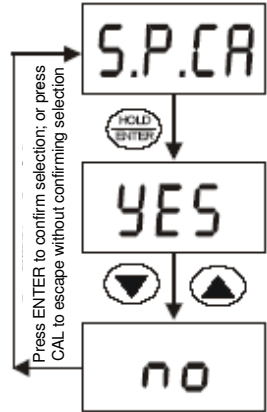


6.8 Restore Factory Default Values

This function allows you to reset all parameters to factory default settings. This clears all calibration data and any other setup functions you might have changed.

IMPORTANT: Once activated, the meter's settings and calibration data will be erased and always exercise caution as meter reset is not reversible.

- 1) Enter the advanced setup as described in Section 6.1.
- 2) Press ▲ or ▼ key until [UrSt] appears on the LCD. Press **ENTER** key.
- 3) Press ▲ or ▼ key to select either [**Yes**] or [**no**].
- 4) Press **ENTER** key to select.
- 5) The meter will go back to measurement mode after the switch-on initialization as shown in figure on Section 3.6.

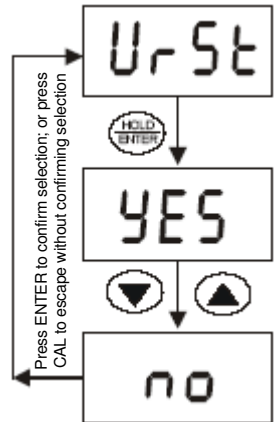

7 PROBE CARE AND MAINTENANCE

Keep the conductivity probe clean. Rinse the probe twice, and gently swirl it while you take readings. For best accuracy, soak a dry probe for at least 5 to 10 minutes or longer before calibration. Rinse the probe with deionized or tap water before storing. Never scratch the bands with a hard substance. Do not strike the probe against any hard surface.


Do not immerse the probe in oily solutions. Clean the electrode thoroughly by stirring it in a mild detergent bath or isopropyl alcohol. Wipe the probe with a soft tissue paper. Rinse thoroughly in tap water and then in deionized water. Recalibrate the meter after cleaning the probe.

The conductivity probe, which is included with your meter, features a removable probe guard to make cleaning easy.

To remove probe guard :-

- 1) Grip yellow probe guard and twist clockwise. The locking notch will release.
- 2) Slide probe guard off end of probe.

NOTE: Remember to re-attach the probe guard prior to taking readings. Failure to do so could result in erroneous readings.

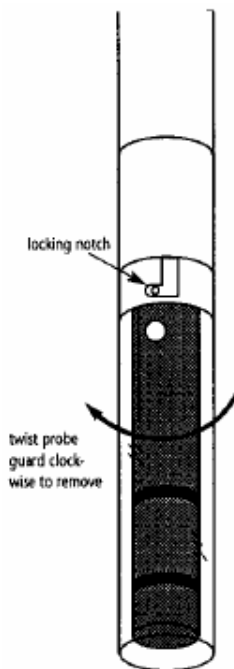

8 TROUBLE-SHOOTING GUIDE

Problem	Cause	Solution
Power On but No Display	<ul style="list-style-type: none"> a) Batteries not in place. b) Batteries not in correct polarity (+ and – position). c) Weak batteries. 	<ul style="list-style-type: none"> a) Check that batteries are in place and making good contact. b) Re-insert batteries with correct polarity. c) Replace batteries.
Unstable Readings	<ul style="list-style-type: none"> a) Air bubbles in probe. b) Dirty probe. c) Probe not deep enough in sample. d) External noise pickup or induction caused by nearby electric motor. e) Broken probe. 	<ul style="list-style-type: none"> a) Tap probe to remove bubbles. b) Clean the probe and re-calibrate. c) Make sure sample entirely covers the probe sensors. d) Move or switch off interfering motor. e) Replace probe.
Slow Response	<ul style="list-style-type: none"> a) Dirty / Oily probe. 	<ul style="list-style-type: none"> a) Clean probe. See 'Probe Care & Maintenance'.


9 ERROR MESSAGES

LCD Display	Indicates	Cause	Solution
'LO' Indicator Appears	Low battery level	Need new batteries or battery connection is bad.	Clean battery contacts. Replace batteries with fresh ones, noting polarity.
Err 1	Conductivity calibration error	Calibration point is outside the $\pm 40\%$ window in the auto-calibration.	Check the value of the conductivity calibration solution. Switch to manual calibration mode and calibrate again. If message persists, return unit.
Err 2	Temperature calibration error	Auto calibration is performed outside the temperature range (0 – 50°C).	Check the temperature and make sure that it is within the acceptable range. If message persists, return unit.
Err 3	Conductivity calibration error	Calibration point is within 10% of the measurement range in the manual calibration mode.	Check the value of the conductivity calibration solution. If message persists, return unit.

If an error message appears, switching off the meter and switching it on again may eliminate the error message. Refer to diagram on right.

If error persists, or the meter shows incorrect values, return the meter.

For a complete diagram of the display see Page 3.

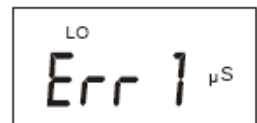

10 SPECIFICATIONS

SPECIFICATIONS	DESCRIPTIONS	PT 142	PT 144
Conductivity Range	0 - 20.00, 200.0, 2000 μ S/cm; 0 - 20.00, 200.0 mS/cm	●	
Resolution	0.01, 0.1, 1 μ S/cm; 0.01, 0.1 mS/cm		●
Accuracy	\pm 1% FS		●
TDS Range	0 - 10.00, 10.0 - 100.0, 100 - 1000 ppm; 1.0 - 10.00, 10.00 - 100.0, up to 200 ppt depending on the TDS factor setting		●
Resolution	0.01, 0.1, 1 ppm; 0.01, 0.1 ppt	●	●
Accuracy	\pm 1% FS	●	●
Temperature Range	-10.0 to 110.0°C	●	●
Resolution/Accuracy	0.1 °C / \pm 0.5 for °C	●	●
Cell Constant	0.1, 1.0, 10.0 (Selectable)	●	●
Temperature Compensation	Automatic / Manual (from 0 - 50°C)	●	●
Temperature Coefficient	0.0 - 3.0% / °C	●	●
Normalization Temperature	20.0°C and 25.0°C (Selectable)	●	●
Conductivity to TDS Conversion factor	0.4 - 1.0		●
Number of Calibration Points	5 (Maximum One per Range)	●	●
Auto and Manual Ranging		●	●
HOLD Function		●	●
Auto Power Off	20 minutes after last key operation	●	●
Inputs	BNC for Conductivity and Phono-Jack for Temperature	●	●
Display	Single Custom LCD	●	●
Power Requirements	4 'AAA' Batteries	●	●
Battery Life	> 100 hours	●	●
Dimension / Weight	14 x 7 x 3.5 cm / 200g	●	●

11 ACCESSORIES

REPLACEMENT PROBES			
PT 142/1	Conductivity / TDS Probe with ATC Sensor		
STANDARD SOLUTIONS FOR CONDUCTIVITY / TDS		CONDUCTIVITY	TDS
PT 142/2	High Range	12.88 mS	8.63 ppt (8630 ppm)
PT 142/3	Mid Range	1412 μ S	946 ppm
PT 142/4	Low Range	74 μ S	50 ppm
