


camGUARD™

Waste Collection & Filter Replacement Service

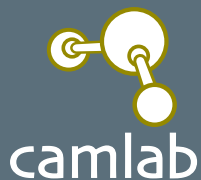

Camlab Ltd
24 Norman Way
Over
Cambridgeshire
CB24 5WE

Contact us:
Tel: 01954 233 110
Fax: 01954 233101
Email: camguard@camlab.co.uk
Web: www.camlab.co.uk

camGUARD™
Hazardous Waste Registration
Ref: NYY846 August 2013
Waste Carbon filter Storage Point:
Camlab Ltd

Safe, environmentally responsible disposal of your used laboratory carbon filter waste.

If hazardous waste is mismanaged it has the potential to cause harm to the environment and human health. As a result, strict controls apply from the point of its production, to its movement, management and disposal.



Combine collection service with waste disposal and a camGUARD™ service contract for the complete solution.

The camGUARD™ system provides a hazardous waste collection and disposal service for saturated carbon filters used in laboratory fume hoods and storage cabinets.

By trapping and binding harmful chemical vapours these filters are key to the protection of laboratory personnel. They need timely replacement and controlled disposal as hazardous waste once saturated.

The camGUARD™ reminder service provides a safe and economical method for the safe replacement and correct disposal of these types of filters. Camlab is authorised to collect, store and consolidate 'used' carbon filters as part of the Non-Waste Framework Directive (NWF) exemptions and as such does not contravene S33(1)a Environmental Protection Act 1990.

This approval by the Environment Agency allows Camlab to operate as a 'Collection Point' for all brands. This economic consolidation of used carbon filters ensures controlled disposal and provides a significant cost saving to the user.

We only use appropriate disposal facilities

Incorrect disposal under the wrong conditions can release harmful chemicals back into the environment.

The camGUARD™ service ensures strict control of disposal to protect the environment. We take great care to ensure complete traceability of the filter and ensure the correct disposal at approved high temperature incineration facilities. Only by complete destruction of the saturated filter at high temperature over 1,100°C can the correct disposal be guaranteed. This is necessary to break down the bound chemicals into non hazardous components. Proof of destruction certificate is available on request.

We handle the paperwork so you don't have to

As a 'Collection point' we take on the responsibility as the waste producer. We provide a fully traceable service producing the hazardous waste consignment notes as well as maintaining records required by the Environment Agency for up to three years.

What should I do when my laboratory filters need replacing?

Register with our camGUARD™ collection service by calling 01954 233110 or emailing camguard@camlab.co.uk. This gives you access to the Asura® range of replacement filters which cover all major brands of hoods and enclosures at special prices.

We can then validate your order and send the filter collection pack documentation for you to complete.



How to get started.....



1 Order your camGUARD™ filter collection pack


Order your replacement filter and the required disposal service.
Please note: for each consignment of filters you will also need to order a single collection note.

Disposal Service

1181278	Contract Customer Asura® Filter Disposal Service including free collection pack	£26.00* per filter
1181277	Non-Contract Customer Asura® Filter Disposal Service including free collection pack	£49.99 per filter
1181279	Duty of care declaration note	

*You can benefit from contract pricing immediately by taking out a new Asura Service Contract. Our engineers can also perform the filter change during the service visit required, as long as the filters have been ordered and delivered to your site in advance.

1168051	Asura® Service Contract Includes service and filter replacement for one hood with a single visit per annum. **Pricing will vary for more than one fume hood, location or where more frequent visits are required. Price does not include filter cost.	from £POA**
1194594	Copy of Proof of Destruction Certificate	£30.00
	FREE camGUARD™ reminder service registration	£FREE

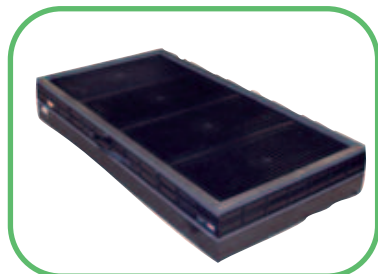

Having ordered your collection service you can also order your replacement Asura carbon filters which will be delivered in reusable packaging.

Please contact us for a quotation.
01954 233110
camguard@camlab.co.uk.

2 Your new filter and camGUARD™ collection pack arrives.

The new filter will arrive in our special packaging which is designed to be opened easily and then re-used for the return filter. If you are using separate packaging supplied by us or need to use your own packaging for any reason please ensure you follow the next steps as closely as possible.




New filter is despatched with a collection notification form and labels for return



3 Pack your used filters in a suitable bag and box

After fitting your new filters, place the used filters in the plastic wrap and seal the bag to ensure it is waterproof using adhesive tape. Then put the bagged filter back into the original carton and seal. If you are fitting your own pre-purchased filters you will need to find a suitable bag and box to return your used filters.

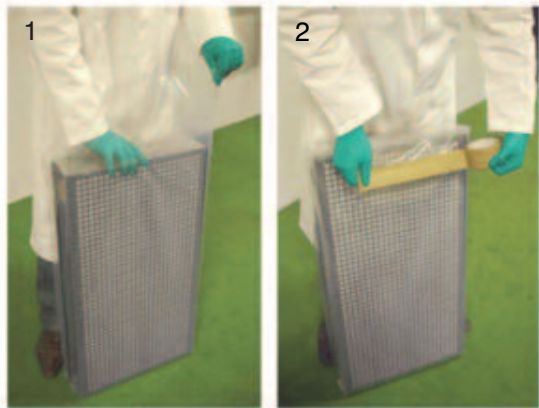

Place the used filter in a plastic bag and seal.



Re-pack in the original box.

4 Label the packages and request collection

Your camGUARD™ pack will come with the correct number of pre-completed transport labels. Simply affix one label to each packed filter and notify Camlab that your filters are ready for collection by faxing the Collection Notification Form to 01954 233 101, or emailing to camguard@camlab.co.uk

NOTE: Please make sure you complete the green sections of the form in full.



Affix label to the carton

camlab camGUARD from Camlab asura FILTERS

COLLECT FROM
Company:
Department of:
Street:
City:
Postcode:

RETURNS NOTE NUMBER
RN 12345
Date: 30/072011

PLEASE FAX THIS FORM TO
01954 233 101

OR EMAIL TO
camguard@camlab.co.uk

INTERNAL USE ONLY
REQUEST RECEIVED: COLLECTION REQUESTED:

5 Paperwork & traceability

As soon as we have received your filters at our consolidation depot we will return the appropriate collection paperwork to you.

A copy will be kept on file for 3 years to provide full environmental traceability of the disposal.

camlab camGUARD from Camlab asura FILTERS

DELIVERY/COLLECTION ADDRESS
Company:
Department:
Street:
City:
POSTCODE:

COLLECTED FROM
Customer:
Department:
Date: 31/07/2011
Returns Note: RN12345

COLLECTION POINT
Camlab Ltd
24 Norman Way Industrial Estate
Over
Cambridgeshire
CB24 5JUE

COLLECTION DECLARATION
I hereby declare that the above information is true and correct to the best of my knowledge and belief. I understand that any false information provided may be a criminal offence under the Environmental Protection Act 1990 and I may be liable to prosecution and/or a fine.
Signature: _____
Name: _____

Denial Ltd
24 Norman Way Industrial Estate, Over, Cambridgeshire, CB24 5JUE
Tel: 01954 233 110 Fax: 01954 233 101 Email: camguard@camlab.co.uk Web: www.camlab.co.uk
Established 1980 Registration Number: 48424