


Nickel 7

Níquel

Test 0-61

Art.-Nr. 985 061

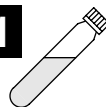

Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.1 - 7.0 mg/l Ni²⁺


470 nm

Method(e) / Método


0611

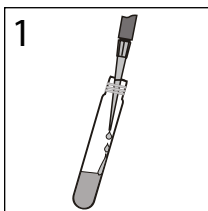


5.0 ml Probe
Sample
Echantillon
Muestra

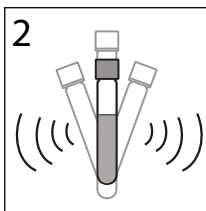


Schütteln
Shake
Agiter
Agitar

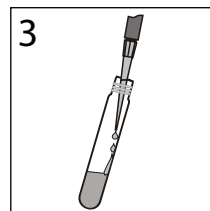


1.0 ml R2

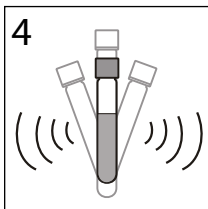


Schütteln
Shake
Agiter
Agitar

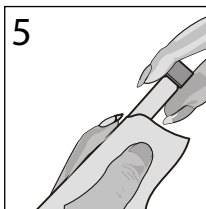

Säubern
Clean
Nettoyer
Limpiar

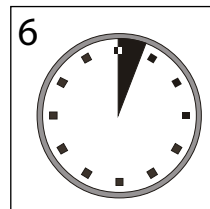

3'00 min

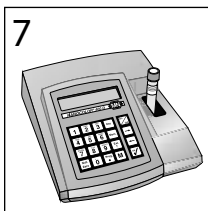

Messung
Measurement
Mesure
Medición