




Built to withstand the harshest working conditions, Aqua washdown scales are rated IP67 and feature a large removable stainless steel weighing pan. Ideal for food production, catering and industrial environments, precision load technology enables improved performance and greater efficiency during weighing and counting tasks. The dual 20mm-high digit front and rear displays enhance visibility and productivity for the user. Checkweighing is simple, with separate high, ok and low coloured LED indicators that change colour to alert the operator if a sample is under, over or within preset limits. Power is supplied by AC adapter or rechargeable battery, and auto power-off conserves battery life.

Get the Aqua wet without worry, the IP67 scale performs effortlessly and accurately no matter how rigorous the application



ABW 32

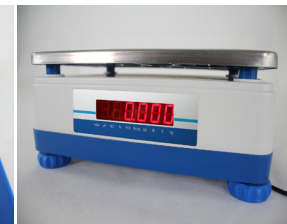

Features

- IP67-rated for protection against water and dust
- External calibration allows for verification and adjustment with weights
- Sealed keypad protects against dirt and spills
- Overload protection helps prevent damage to internal components
- Solid ABS housing is durable for industrial use
- Heavy-duty rubber feet ensure stability of unit
- Bright LED display for easy viewing in most environments
- Checkweighing LEDs clearly show under, over and within acceptable limits
- Front and rear 20mm LED displays allow for dual use
- Zero-tracking feature ensures display returns to zero reading
- The auto power-off feature conserves battery life
- Rechargeable battery included for operation almost anywhere
- Stabilization time of within 1 second
- AC adapter included

Applications

- Weighing
- Parts counting
- Percentage weighing
- Checkweighing

About ADAM:


For nearly 50 years, professionals worldwide have relied on Adam Equipment for an extensive selection of dependable and affordable weighing equipment. Adam provides the right balance of speed, performance and value.

Product Details

Model	ABW 4	ABW 8	ABW 16	ABW 32
Capacity	4kg	8kg	16kg	32kg
Readability	0.5g	1g	2g	5g
Repeatability	1g	2g	5g	10g
Linearity	1g	2g	5g	10g
Pan Size	245x180mm			300x210mm
Weighing Units	g, kg, lb, oz, lb:oz			
Stabilization Time (sec.)	1			
Calibration	External			
Display	Dual-display LED with 20mm-high digits			
Power Supply	12VDC 800mA adapter, Internal rechargeable battery			
Operating Temperature	0° to 40°C			
Housing	ABS plastic			
Overall Dimensions	255x290x105mm			300x310x105mm
Net Weight	2.5kg			

Accessories

Item number	Description
3022014061	In-use cover
2020014062	In-use cover (pack of 5)
2020014063	In-use cover (pack of 10)
3002014371	Hard carrying case
303149759	Vegetable Scoop (complete with fitting to scales)
303149760	Confectionery Scoop (complete with fitting to scales)
700660290	Calibration certificate

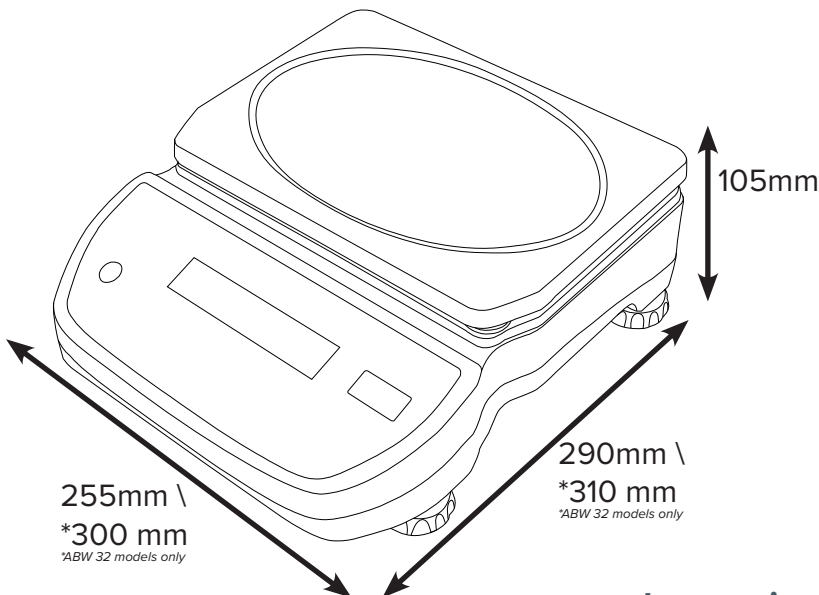

Key Accessories



Scoops

Metal scoops are designed to make portioning and food measurement easy. Two sizes available: small and large.



Hard Carrying Case

Allows the Aqua to be placed inside to protect it during transportation. The hard shell case is durable and comes with a lock so that you can secure the scale in the case as needed. Die cut foam insert provides room for the scale, power supply and top pan.

Supplied by:



24 Norman Way Industrial Estate, Over, Cambridge, CB24 5WE
Tel: +44(0)1954 233 100 Fax: +44(0)1954 233 101
Email: sales@camlab.co.uk Web: www.camlab.co.uk

www.adamequipment.co.uk

Adam Equipment follows a policy of continuous improvement and reserves the right to change specification and pricing without notice. Adam Equipment's standard terms and conditions of sale apply. Errors and omissions excepted.