


Confident Cooling™

A New Way To Chill

Grant

"Grant Thermocirculators are widely used in the Chemistry Dept at Imperial College, London".

Prof. T Welton - Imperial College, London

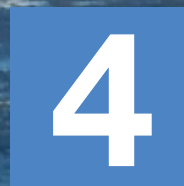

Confident Cooling™ Guarantee

4 Year Warranty

LTC Refrigerated / Heating circulating baths:



Confident Cooling™ for most applications, including:

Rotary evaporation

Reaction vessel cooling

Condensers

Electrophoresis



Excellent cooling power: 340W and 900W options



Ready to use kits



Dispatched within 5 working days



4 year best in class warranty



Manufactured in Cambridge, UK



Expert support by our dedicated Confident Cooling™ Support Team



Full refund if not entirely satisfied



Confident Cooling™ Support Team helpline:

+44 (0) 1763 264 778

confident.cooling@grantinstruments.com

www.grantinstruments.com

Product Specifications

	LTC2 General Purpose Optima™	LTC4 Advanced Optima™
		
Ready to go kits	Thermostat, refrigerator, insulated tubing and clips provided to form a ready-to-use system	
Warranty	4 year warranty with no product registration needed	
Dispatch	Within 5 working days from date of order	
Pricing	Highly competitive	
Manufacture	United Kingdom	
Number of pre-set temperatures	3	
Programmable	-	1 Program / 30 segments via Labwise™
Cooling power @ 20°C*	340W	900W
Display	Simple dial and button user interface	Full colour TFT screen with remote control option in 5 languages
Temperature range	-25°C to 100°C	-30°C to 100°C
Temperature alarms	User defined high temperature alarm	User defined high and low temperature alarms
User calibration	2 point user calibration	5 point user calibration

Product Features



*Confident Cooling™
Guarantee*



Ready to go kits



*4 year warranty with no
registration needed*



*High cooling power for
multiple connections**



*Tubing can be cut to
specified length easily*



*A choice of tubing pump
connectors included as standard*



*High pressure pump option can be
fitted for demanding applications*

Confident Cooling™ Support Team helpline:

+44 (0) 1763 264 778

confident.cooling@grantinstruments.com

* Typical values

+ Shown with Asynt ReactoMate

www.grantinstruments.com