



Well within Reach

The Eppendorf Pipette Holder System



Pipette Stand 2 and Charger Stand 2



Protect your precious liquid handling device from possible damage, contamination or chemicals spilled on the bench. The stands for single instruments have a robust design and securely hold one manual or electronic device. With the Charger Stand 2, your Eppendorf Xplorer®, Eppendorf Xplorer plus or Multipette® E3/E3x is always fully charged and ready for use. The completely new Pipette Stand 2 carries one Eppendorf Multipette M4. This stand is also perfect for one single or multi-channel Eppendorf Research® plus or Eppendorf Reference® 2 pipette when used with the corresponding pipette holder.



Perfection Down to The Smallest Detail

Rotatable holders carry all current manual Eppendorf pipettes and most predecessors.



Pictograms clearly label pipette holders and charger shells for easy assignment in the lab.



Fast interchangeable pipette holders and charger shells, no tool required.



Downwards compatible holders for Eppendorf Research and Reference manual pipettes.



Large carrying handle for easy transport of carousels from bench to bench.



Large rubber feet protect carousels and stands from liquids spilled on lab bench.



Pipette Holder 2 and Charger Shell 2

Choose from a broad range of holders for all current manual and electronic Eppendorf liquid handling instruments and most predecessors. Place your pipette and Multipipette where you need them: Mounted to a wall, on the shelf above your bench or inside

a biological safety cabinet. In addition, all holders can be inserted into Pipette Carousel 2. Charger Shells for installation into Charger Carousel 2 are available as well.

Robust electrical contacts designed for long operational lifetime without the need for replacement.




V-shaped charging contacts ensure stable electrical connection between charger shell and electronic device.



Cord wrap functionality for unused power cord keeps your lab bench nice and tidy.



A small footprint with a diameter of 18.5 cm makes stands and carousels ideal when space is at a premium.

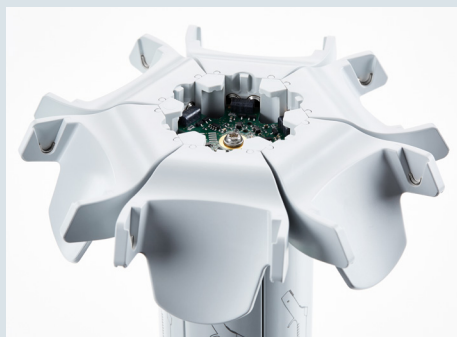

Intelligent charging electronics keeps an eye on up to six Eppendorf Xplorer pipettes or electronic Multipipettes.



Quick-release magnetic power connector protects pipettes and Multipipettes from falling.

Pipette Carousel 2 and Charger Carousel 2

Keep all your Eppendorf pipettes and Multipettes well within reach. The Pipette Carousel 2 for manual pipettes and the Charger Carousel 2 for electronic devices offer more flexibility and more capacity than ever before.

The Eppendorf Pipette Carousel 2 carries up to six manual single or multi-channel Eppendorf pipettes. Single pipette holders can be mounted to the carousel in two orientations turned by 180 degrees to each other to store all current and most predecessor manual pipettes. This includes Eppendorf Research, Research plus, Reference, Reference 2 and Biomaster®.

The Charger Carousel 2 offers 50% more capacity than before and holds up to six electronic Eppendorf Xplorer or Xplorer plus pipettes. When used with optional charger shells the carousel can carry up to six Eppendorf Multipette E3/E3x or Multipette stream/Xstream. In addition to electronic devices, the Charger Carousel 2 is also perfect for manual Eppendorf pipettes when used with optional pipette holders.

More flexibility

- > Carousels carry all kinds of Eppendorf pipettes and Multipettes
- > All manual Eppendorf pipettes fit into one rotatable holder



More capacity

- > Charger Carousel 2 for up to six electronic devices
- > Carousels carry both manual and electronic devices

More convenience

- > Interchangeable pipette holders, no tools required
- > Cord wrap functionality for unused power cord

More user safety

- > Magnetic power connector for fast and easy connection
- > Large rubber feet protect carousels from liquids spilled on bench top



> Additional information available at:
www.eppendorf.com/pipetteholder


Ordering information

Description	Order No.
Pipette Carousel 2 , for 6 Eppendorf Research®, Eppendorf Research® plus, Eppendorf Reference®, Eppendorf Reference® 2 or Biomaster®, additional pipette holders are optionally available	3116 000.015
Charger Carousel 2 , for 6 Eppendorf Xplorer® or Eppendorf Xplorer® plus, mains/power adapter included, additional charger shells and pipette holders are optionally available	3116 000.023
Charger Stand 2 , for one Eppendorf Xplorer® or Eppendorf Xplorer® plus, operated with mains/power adapter supplied with Eppendorf Xplorer® or Eppendorf Xplorer® plus	3116 000.031
Charger Stand 2 , for one Eppendorf Multipette® E3/E3x or Multipette® stream/Xstream, operated with mains/power adapter supplied with Eppendorf Multipette® E3/E3x or Multipette® stream/Xstream	3116 000.040
Pipette Stand 2 , for one Eppendorf Multipette® M4, without charging functionality, additional pipette holders are optionally available	3116 000.058
Pipette Holder 2 , for one Eppendorf Research®, Eppendorf Research® plus, Eppendorf Reference®, Eppendorf Reference® 2 or Biomaster®, for Pipette Carousel 2 and Charger Carousel 2 or wall mounting, sticky tape included	3116 000.112
Pipette Holder 2 , for one Eppendorf Xplorer® or Eppendorf Xplorer® plus, for Pipette Carousel 2 or wall mounting, sticky tape included, without charging functionality	3116 000.120
Pipette Holder 2 , for one Eppendorf Multipette® E3/E3x or Multipette® stream/Xstream, for Pipette Carousel 2 or wall mounting, sticky tape included, without charging functionality	3116 000.139
Pipette Holder 2 , for one Eppendorf Multipette® M4, for Pipette Carousel 2 and Charger Carousel 2 or wall mounting, sticky tape included, without charging functionality	3116 000.147
Charger Shell 2 , for one Eppendorf Xplorer® or Eppendorf Xplorer® plus, for Charger Carousel 2, with charging functionality	3116 602.007
Charger Shell 2 , for one Eppendorf Multipette® E3/E3x or Multipette® stream/Xstream, for Charger Carousel 2, with charging functionality	3116 603.003

Learn more about the new Pipette Holder System

The Eppendorf Pipette Holder System provides highly flexible carousels, stands and wall mount devices for all current manual and electronic Eppendorf liquid handling instruments and most predecessors.

> www.eppendorf.com/pipetteholder

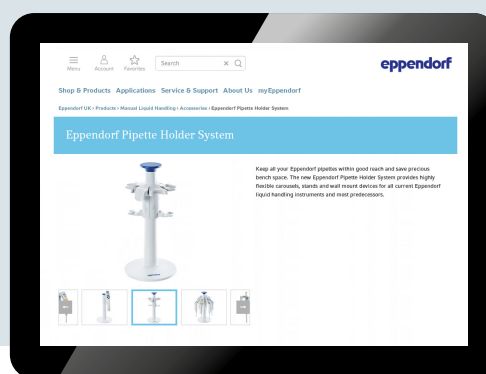

24 Norman Way Industrial Estate, Over, Camb, CB24 5WE

Tel: +44(0)1954 233100 Fax: +44(0)1954 233101
Email: sales@camlab.co.uk Web: www.camlab.co.uk