


RANGE 2**FM - NFPA - OSHA**

FM SAFETY CABINETS for flammable products

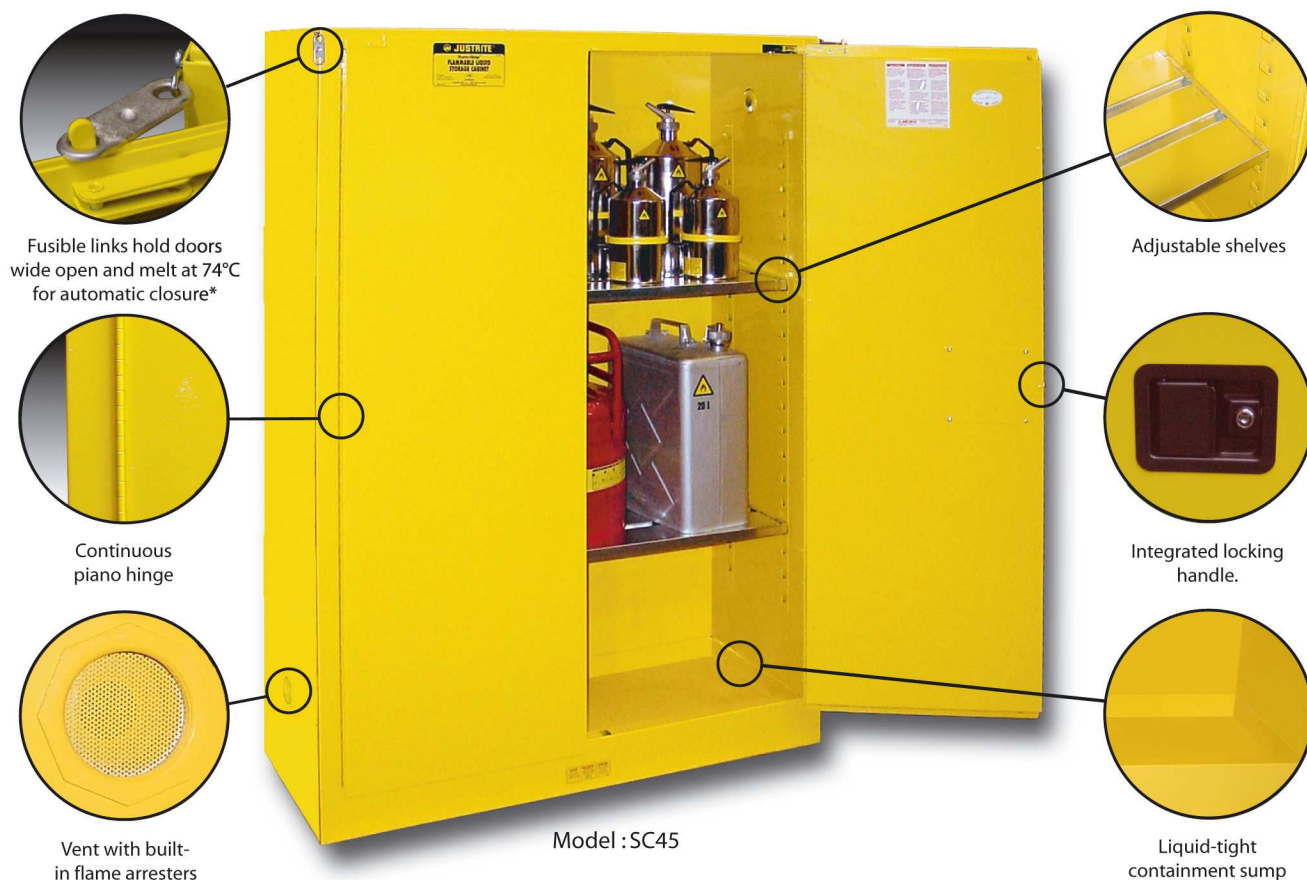

Model : SC45

DIRECTIVES & NORMS

- In accordance with NFPA (National Fire Protection Association), OSHA (Occupational Safety and Health Administration), and FM (Factory Mutual international) regulations.

ADVANTAGES

- Sturdy construction.
- Integrated locking handle.
- Cabinets with self-closing doors. The fusible links hold doors wide open and melt at 74°C for automatic closure*
- Continuous piano hinge provides smooth closure.
- Many models are available.

All FM approved cabinets are recommended by the U.S. insurance companies.

PASSIVE SECURITY

- Built in double wall steel with 38 mm of insulating air space.
- Yellow epoxy paint.
- Adjustable shelves.
- 50 mm liquid-tight containment sump.
- Built in grounding connector.

ACTIVE SECURITY

- Lateral high and low vent Ø 2" (50mm) with flame arresters.
- Continuous piano hinge.
- Manual or self-closing doors.
- 3 points lock with fail-safe closing mechanism and locking handle.
- Key locking system.

OPTIONS

- Additional shelf.
- Absorbing carpets (ref. JF10).
- Ventilation or filtration boxes.

* On self-closing models.

RANGE 2

FM - NFPA - OSHA



Ref.	Door types	Description	External dimension H x W x D (mm)	Internal dimension H x W x D (mm)	Storage capacity (L)	Retention capacity (L)	Number of shelves	W. (kg)	
MB04	Manual	Under-bench cabinet with 1 door	610 x 430 x 430	465 x 345 x 360	15	6	1	30	
SB04	Self-closing							32	
MB12	Manual	Working cover cabinet with 1 door	940 x 610 x 460	795 x 510 x 390	45	10	1	65	
SB12	Self-closing							67	
ES12	-	Additional shelf for 45L cabinet							3
MB18	Manual	Under-bench cabinet with 2 doors	610 x 1090 x 460	515 x 1005 x 390	64	20	1	75	
SB18	Self-closing							78	
ES18	-	Additional shelf for 64L cabinet							4
MB22	Manual	Tall cabinet with 1 door	1650 x 585 x 460	1560 x 500 x 390	83	10	3	90	
SB22	Self-closing							92	
ES22	-	Additional shelf for 83L cabinet							2
MB30	Manual	Working cover cabinet with 2 doors	1120 x 1090 x 460	1025 x 1005 x 390	113	20	1	105	
SB30	Self-closing							108	
ES30	-	Additional shelf for 113L cabinet							4
MC45	Manual	Tall cabinet with 2 doors	1650 x 1090 x 460	1560 x 1005 x 390	170	20	2	140	
SC45	Self-closing							143	
ES45	-	Additional shelf for 170L cabinet							4
MC60	Manual	Tall cabinet with 2 doors	1650 x 860 x 860	1560 x 770 x 780	228	30	2	170	
SC60	Self-closing							173	
ES60	-	Additional shelf for 228L cabinet							8
MC90**	Manual	Tall cabinet with 2 doors	1650 x 1100 x 860	1560 x 1020 x 780	340	40	2	220	
SC90**	Self-closing							223	
ES90**	-	Additional shelf for 340L cabinet							10

** Model not on stock and not shown on the picture.



Some models of Range 2 FM safety cabinets are designed to stack on Range 4 cabinets. Contact-us for more information.

To protect the users health and to improve the life cycle and the safety of your cabinet, we recommend to install an exhausting system. Contact us for additional information.

