



Captair Smart

Ductless filtering fume hoods

Life in the lab becomes Simpler *Safer*



Captair Smart

Ductless Filtering Fume Hoods

Erlab's innovative technology keeps you safe

Erlab's unique Flex Filtration can manage most lab applications

SMART TECHNOLOGY

Simple, intuitive communication by light shows the status of the hood

Sensors to detect filter saturation

Low power consumption. Operates with a standard plug using only 60W


Face velocity is automatically monitored

The 6° slanted sash offers optimal visibility and an ergonomic work space

Large openings for easy access to your work



Download the eGuard App for remote control and monitoring

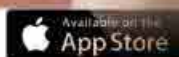

Life in the laboratory becomes simpler and safer

Simpler to use *Safer* to operate

Erlab's 45 years of expertise in research, design and manufacturing of filtering fume hoods guarantee superior filtration that will keep you protected during your chemical processes.

Smart Technology: The new Captair Smart ductless filtering fume hoods are designed with a simple and innovative way of communication. Smart Technology uses simple light to show that the hood is operating safely so you can focus your attention on what is most important: your work.



Erlab's expertise
in filtration



Guided by
the light



Remote access to hood
status, settings and alerts

Simpler to use

- Smart Technology enables visual light communication to indicate the status of the hood. When the light is still the hood is operating safely. To alert you, the light will pulse. A simple light pattern will indicate if :
 - Air flow is compromised
 - Fan Failure
 - Filter breakthrough
- The eGuard App gives you remote control to monitor the hood, change the settings, and delivers safety alerts immediately to your mobile, tablet or PC device.





Safer to operate

- The configurable filtering system above the enclosure can be adapted to filter gases, solvents, powders and particulates keeping the user and the lab protected by safely and efficiently recirculating air within the room and releasing purified air back into the room free of any toxic hazardous vapors or odors
- Meets AFNOR NFX 15 211/ANSI Z9.5-2012 filtration efficiency standard (classes 1 and 2).
- New sensors are added for detection of solvents or acids or formaldehyde.
- New electronic anemometer monitors face velocity
- ESP Program: Usage certification delivered with each unit after approval of your handlings or process by our in-house lab.

321

391

481

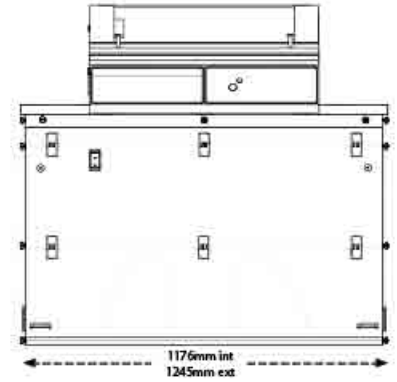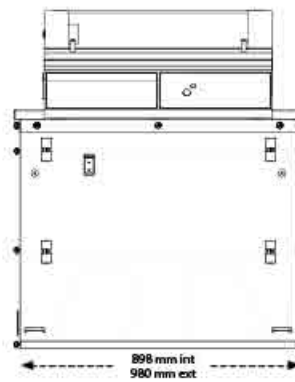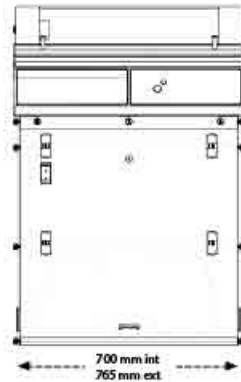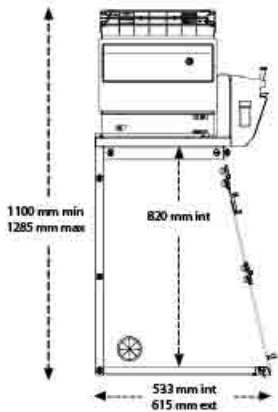

Safety Standards	AFNOR NF X 15-211:2009 : France - BS 7989 : England DIN 12 927 : Germany - EN 1822 : 1998 (HEPA H14 & ULPA U17 Filters) - CE Marking
Air Flow	200 to 230 m ³ /h - 117 CFM to 135 CFM
Air Face Velocity	0,4 to 0,6 m/s - 79 fpm to 118 fpm
Voltage/Frequency	80 - 220 V / 50-60 Hz
Power consumption	60 W
Sash openings	Oblong
Structure	Corrosion resistant electro-galvanized steel coated with anti-acid polymer
Side and front panels	Chemical resistant acrylic
Filtration module	Polypropylene

Features

Smart Technology	Simple communication by LED pulsation system: fan settings, usage timer, fan failure, face velocity, automatic filter saturation detection
Flex Technology	1 column that can be configured to handle liquids, powders, or both
Carbon filtration for gases and vapors	AS : For organic vapors - BE+ : For organic vapors and acid vapors F : For formaldehyde vapors - K : For ammonia vapors
Particulate filtration for powders	HEPA H14 : filtration efficiency: 99,995 % according to MPPS method, EN1822 standard ULPA U17 filtration efficiency: 99,999995 % according to MPPS method, EN1822 standard
eGuard	App for remote control to monitor the hood, change the settings, and deliver safety alerts immediately to your devices (mobile, tablet and PC).
Molecode	Detection sensor for solvents or acids or formaldehydes
Internal lighting	LED lighting > 650 Lux
Electronic Anemometer	Indicates the face velocity of the unit
Chemical Listing	List of approved chemicals

Accessories

Work Surfaces	Trespa® Top Lab ^{plus} , Glass or 304L Stainless Steel
Benches	Rolling or Fixed
Bench equipment	Technical gases outlets, water outlets, front control valves, sink, power sockets (Only compatible with Trespa® Top Lab ^{plus} worktop and fixed bench)
Particulate Pre-filter	40 % filtration efficiency for particles > 0,3 µm
Wasteport	Double-bag with protected housing
Transparent Back Panel	Clear acrylic panel for easy viewing

France
+33 (0) 2 32 09 55 80 | ventes@erlab.net

USA
+1 800-964-4434 | captairsales@erlab.com

China
+86 (0) 512 5781 4085 | sales.china@erlab.com.cn

Malaysia
+60 (0) 7 3 555 724 | erlab@tm.net.my

Germany
0800 330 47 31 | verkauf@erlab.net

United-Kingdom
+44 (0) 1722 341 940 | salesuk@erlab.net

Italy
+39 (0) 2 89 00 771 | vendite@erlab.net

Spain
+34 93 673 24 74 | ventas@erlab.net

