

THERMOSTAT DRB200 – THE SIMPLE SOLUTION FOR ALL DIGESTIONS



PRODUCT INFORMATION

- Laboratory analysis
- Dry thermostat
- DRB200
- **For 16 mm cuvettes**

For standard and special digestions

A digestion step often has to be carried out before key parameters can be determined. For these parameters, the DRB200 thermostat ideally complements the Hach reagents (COD, TOC, Test'N'Tube). Two separately controllable heating blocks enable cuvettes and reaction vessels to be digested at identical or different temperatures and time settings.

The DRB200 has a digital timer with an automatic switch-off and acoustic signal. Two transparent splash protection lids close the thermostat while it is heating. The integrated anti-overheating feature and the insulated external shell provide additional safety.

Great flexibility

Pre-programmed for all standard digestions and freely programmable for user-specific digestions

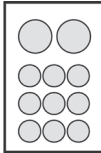
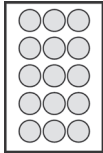
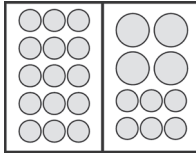
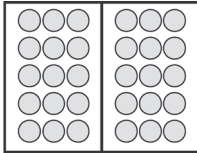
Simple handling

Large, easily readable display and one-key operation for standard digestions

Excellent reproducibility

Temperature stability

Technical Data

Type	DRB200-1	DRB200-1	DRB200-2	DRB200-2
Order no.	LTV082.53.30001 (115 VAC) LTV082.52.30001 (230 VAC)	LTV082.53.40001 (115 VAC) LTV082.52.40001 (230 VAC)	LTV082.53.42001 (115 VAC) LTV082.52.42001 (230 VAC)	LTV082.53.44001 (115 VAC) LTV082.52.44001 (230 VAC)
Heating programs	Pre-programmed for 100 °C, 105 °C, 150 °C, 165 °C and freely selectable 37–150 °C, 1–480 min.			
Heating rate	From 20 °C to 150 °C in 10 min.			
Temperature stability	±1 °C in conformity with EN, ISO, EPA methods			
Power input	115 V / 300 VA 230 V / 450 VA		115 V / 600 VA 230 V / 900 VA	
Weight	2 kg (4.4 lbs)		2.8 kg (6.2 lbs)	
Dimensions	250 × 145 × 310 mm (W×H×D)			
Schematic representation				
Number of cuvettes	9 × 16 mm diameter 2 × 20 mm diameter	15 × 16 mm diameter	21 × 16 mm diameter 4 × 20 mm diameter	30 × 16 mm diameter



DRB200-1 for 9 cuvettes (16 mm ø) and 2 reaction vessels



The DRB200-2 features two separately controllable heating blocks

Applications	Temp. [°C]	Time [min.]
COD	150	120
Total nitrogen	100 / 105*	60 / 30
Total phosphorus	100 / 105*	60 / 30
Metals (lead, cadmium, copper, iron, nickel, zinc)	100	60
Total chromium	100	60
TOC	105	120
Trihalogenmethane	100	8
User-specific programs	37–150	1–480

* Digestion temperature and time depend on the used reagents. Subject to change without notice.

HACH COMPANY World Headquarters: Loveland, Colorado USA

United States: 800-227-4224 tel 970-669-2932 fax orders@hach.com
 Outside United States: 970-669-3050 tel 970-461-3939 fax int@hach.com
hach.com

Printed in U.S.A.

©Hach Company, 2017. All rights reserved.

In the interest of improving and updating its equipment, Hach Company reserves the right to alter specifications to equipment at any time.



24 Norman Way Industrial Estate, Over, Cambridge, CB24 5WE
 Tel: +44(0)1954 233 100 Fax: +44(0)1954 233 101
 Email: sales@camlab.co.uk Web: www.camlab.co.uk

DOC062.53.00362.May17



Be Right™