




The **installation drawing** is available under [www.kirsch-medical.com](http://www.kirsch-medical.com) for each model in the download area.

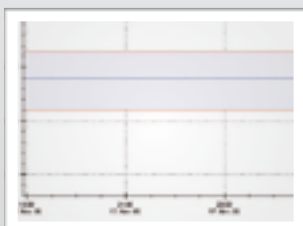

**Constant temperature**


Maintains the interior temperature (maximum: -22 °C) even in the case of greater ambient temperatures of up to +32 °C.



**Suitable for integration**

Compact freezer designed for the integration into furniture. Can also be used as a standalone device.

# FROSTER- LABEX<sup>®</sup>-70 II 3 G T6

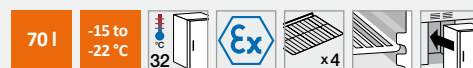

- **External housing** made from galvanised sheet steel (rust-proof), with white, anti-scratch powder coating. Length of the plug cable: approx. 1.8 m.
- **Interior** made from impact-resistant, white plastic with moulded-in shelves.
- **Interior consists** of four plastic-coated shelves.
- **Extra-thick 55-mm insulation**, made from high-quality, compression-moulded and environmentally friendly material. Energy saving.
- **Lockable door** with easy-to-replace plastic magnetic seal frame.
- **Door hinge** on the right-hand side (by default, see illustration), or the left-hand side, can be retrofitted.
- **Manual defrosting.**
- **Mechanical controller** enables variable temperature control. The temperature is maintained automatically, irrespective of changing external temperatures, provided that the ambient temperature is at least approx. 3 °C above the set interior temperature.
- **Statically ventilated refrigerating machine**, hermetically sealed, energy saving, low-noise, easy to service.
- **Door frame heater** prevents icing of the door.

## Specification

Capacity	70 litres
Temperature setting	approx. -15 °C to -22 °C
Voltage	220 – 240 V, 50/60 Hz
Power consumption	100 watts
Normal consumption	1.6 kWh/24 h
Admissible ambient temperature	from +10 °C to +32 °C
Heat emission (max.)	129 watts
Exterior dimensions (including distance from wall)	w x d x h = 47 x 51 x 74 cm
Interior dimensions	w x d x h = 35 x 41 x 55 cm (usable depth above: 9 cm, below: 20 cm less)
Exterior dimensions with door open at 90°	w x d = 47 x 96 cm
Shelfsize	w x d = 35 x 31 cm (upper), w x d = 35 x 20 cm (lower)
Max. load shelf	25 kg
Weight	Net weight 34 kg, gross weight 42 kg

## Optional equipment

- **Decorative door frame** for attaching the customer's decorative plates.
- **Door fitting** (two pieces required).
- **Equipped with castors.**
- **Refrigerating machine fan:** required when the system is used as a built-in freezer if static ventilation is insufficient.
- **Temperature documentation:** see detailed information on page 56.