


MS-H-Pro+ LCD Digital Magnetic Hotplate Stirrer

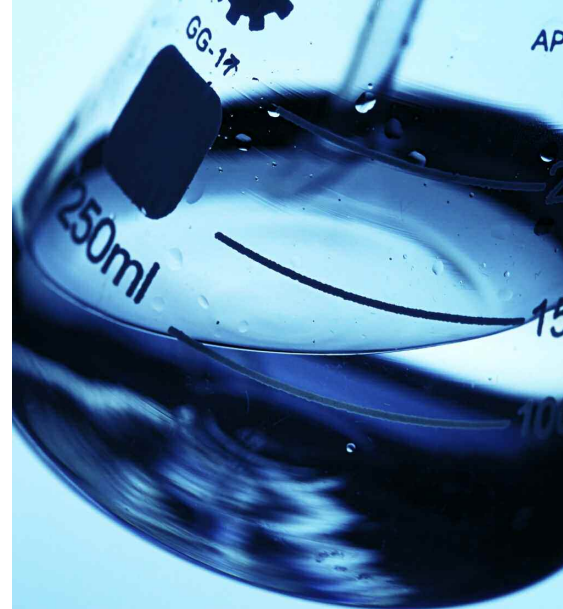

Features

- Stainless steel or ceramic plate
- DC brushless motor
- Heating Range RT-340°C
- LCD display
- HOT warning indicator

The Camlab Choice MS-H-Pro+ is widely used in chemical synthesis, bio-pharmaceutical and clinical laboratories. PID temperature technology enables the unit to rapidly reach optimum operating temperature. The stainless steel hotplate is impact and corrosion resistant with good uniformity of heating.


The MS-H-Pro+ can be fitted with the PT1000 temperature sensor (as illustrated) which is automatically identified when connected for real time control.

The DC Brushless motor is maintenance free and very quiet in operation.

The M-H-Pro+ is available with a ceramic coated or stainless steel hotplate.

Specifications

Part Number	1189931
Frequency	50/60Hz
Max stirring qty	20L
Speed range	100-1500 rpm
Heat output	500W
Safety temp	360°C
Dimensions	160x280x85
Weight	2.8Kg
Workplate dimension	135mm

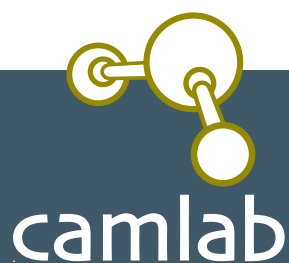



camlab

Camlab Limited, Camlab House, Norman Way, Over, Cambridge, CB24 5WE
Tel: 01954 233110 • Fax: 01954 233101
email: mailbox@camlab.co.uk • www.camlab.co.uk