





















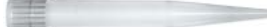


## Socorex® compatibility chart








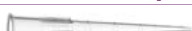














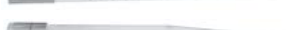
Chart indicates optimal combination of Qualitix® and Socorex® pipettes.

		Acura® manual and Acura® electro												Calibra® digital																			
		815					825 / 826 / 926					835 / 936		855 / 956			822		832		852												
		Up to 10 µL	10 µL / 15 / 20 µL	25 to 100 µL	120 to 200 µL	250 to 1000 µL	0.1 - 2 µL	0.5 - 10 µL	1 - 10 µL (Y)	2 - 20 µL	5 - 50 µL	10 - 100 µL	20 - 200 µL	100 - 1000 µL	0.2 - 2 mL	0.5 - 5 mL	1 - 10 mL	0.5 - 10 µL	5 - 50 µL	10 - 100 µL	20 - 200 µL	40 - 350 µL	2 - 20 µL	10 - 100 µL	20 - 200 µL	100 - 1000 µL	0.2 - 2 mL	1 - 10 mL	1 - 10 µL	10 - 100 µL	20 - 200 µL		
<b>10 µL Microtips</b>	<b>Tip ref.</b>																																
	3090	●					●	●										●															
	3020	●					●	●																									
	302	● <sup>2</sup>							● <sup>2</sup>																								
	3090	●					●	●																									
	3020	●					●	●																									
<b>20 - 100 µL Microtips</b>	<b>Tip ref.</b>																																
	3080		●						●	●	● <sup>1</sup>																						
	3080			●							●	●	● <sup>1</sup>																				
<b>200 µL Microtips</b>	<b>Tip ref.</b>																																
	3070	●	●	●					●	●	●	●	●							●	●	●											
	3290	●	●	●					●	●	●	●	●							●	●	●											
	328																						●	●	●					●	●		
	3080				●					●	●	●	●																				
	3070		●	●									●																				
	3080			●									●																				
<b>300 - 350 µL Microtips</b>	<b>Tip ref.</b>																																
	308			●								●	●										●	●						●			
	308			●									●											●						●			
<b>1000 µL Microtips</b>	<b>Tip ref.</b>																																
	3070				●								●																				
	3190				●								●																				
	3090				●								●																				
	319				●								●													●							
	3070				●								●																				
<b>2, 5, 10 mL Macrotips</b>	<b>Tip ref.</b>																																
	312														●												●						
	312															●												●					
	312																●												●				

<sup>1</sup> Up to max. tip volume    <sup>2</sup> For volumes > 2 µL

## Other brands compatibility chart

The Qualitix® line also fits major micropipette brands available in the market.

	Brand	Eppendorf					Gilson									
		Such as models: Transferrpette S, Electro					Such as models: Research Plus / Pro					Such as models: Pipetman, Classic, G, L				
		2.5, 10 µL	20, 50, 100, 200 µL	250, 300 µL	1000 (1200) µL	2, 5, 10 mL	2, 10 µL	20, 50, 100, 200 µL	300 µL	1000 µL	5, 10 mL	2, 10 µL	20, 50, 100, 200 µL	250, 300 µL	1000 µL	2, 5, 10 mL
<b>10 µL Microtips</b>	<b>Tip ref.</b>															
	3090	●					●					●				
	3020	●					●					●				
	302	●					●					●				
	3090	●					●					●				
	3020	●					●					●				
<b>20 - 100 µL Microtips</b>	<b>Tip ref.</b>															
	3080						● <sup>1</sup>					● <sup>1</sup>				
	3080	● <sup>1</sup>					● <sup>1</sup>					● <sup>1</sup>				
<b>200 µL Microtips</b>	<b>Tip ref.</b>															
	3070		●					●					●			
	3290		●					●					●			
	328															
	3080		●					●					●			
	3070		●					●					●			
	3080		●					●					●			
<b>300 - 350 µL Microtips</b>	<b>Tip ref.</b>															
	308			●									●			
	308			●												
<b>1000 µL Microtips</b>	<b>Tip ref.</b>															
	3070				●				●					●		
	3190				●				●					●		
	3090				●				●					●		
	319				●				●					●		
	3070				●				●					●		
<b>2, 5, 10 mL Macrotips</b>	<b>Tip ref.</b>															
	312									●						
	312															●
	312															

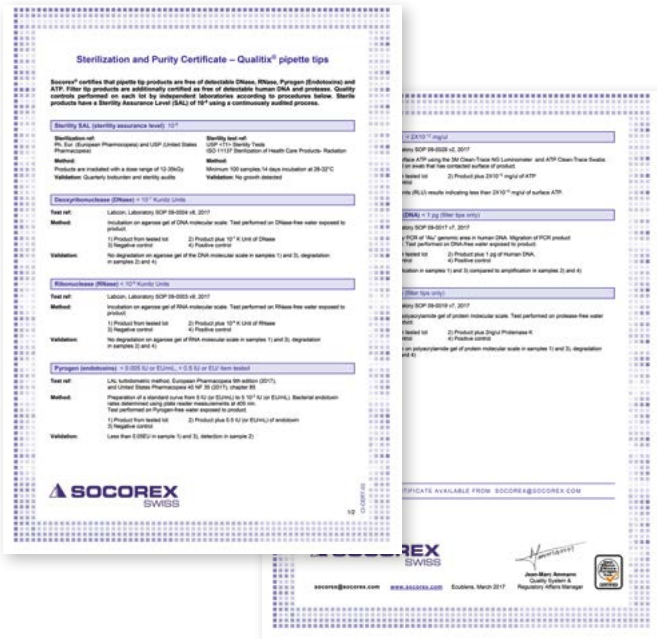
<sup>1</sup> Up to maximal tip volume

Brand names belong to their respective owners.



## Superior quality material

Socorex is aiming at providing reliable pipette tips to ensure consistent analytical results. High purity grade and lot-by-lot consistency guaranteed. Product certificate supplied in each pack documents test methods and detection level. Lot specific certificate available from Socorex and its local partners.



View certificate



### ► First class material

Only virgin, medical grade polypropylene used in tip manufacturing process. Material guaranteed metal-free and non cytotoxic, eliminating possible disturbance.

### ► bioproof level

The Qualitix® line is free of detectable DNase, RNase, Pyrogen (Endotoxins) and ATP. Filter tips are additionally certified free of detectable human DNA and protease.



## Purity benefits

Feature	Advantage	Regular Qualitix®	Low binding tips	Filter tips
Nuclease-free (DNase and RNase)	Safe handling of samples containing genetic material	•	•	•
Human DNA-free	bioproof purity level for any critical pipetting (PCR, qPCR)			•
Pyrogen-free	Recommended in cell culture laboratories	•	•	•
ATP-free (Adenosin Triphosphate)	Beneficial in accurate measurement of cells or bacteria activity	•	•	•
Heavy metal-free	Adequate for trace analysis	•	•	•
Protease-free	Protection of valuable samples containing peptides and proteins		•	•
Sterility	Preserve sample integrity	• <sup>1</sup>	• <sup>1</sup>	•

<sup>1</sup> Sterile items only

### ► Stringent quality controls

Each lot is tested by independent laboratories according to recognized procedures. Sterile products have a Sterility Assurance Level (SAL) of 10<sup>-6</sup> using a continuously audited process.

### ► Standards and compliance

Qualitix® pipette tips conform to the CE directive IVD 98/79 EEC (in vitro diagnostic medical devices).



# Qualitix®

reliable consumables – by excellence

# premium micropipette tips



### All-purpose reagent reservoirs

Well suited for single and multichannel pipetting. V-shape bottom enables aspiration until last reagent drops. Made of PP, autoclavable, re-usable.



More details



### 8- and 12-channel reservoirs

Made of PS material, provide ideal solutions for multichannel pipette applications. Supplied in bulk pack (non-sterile) or single wrapped (sterile).



### Straw tips for 1:10 dilution pipette

Fitting the Acura® 810 model, the long and thin shape straw tips allow pick-up of particle containing liquids in Stomacher® bags and deep, narrow reservoirs. Supplied sterilized in re-sealable bags.



More details



### Glass Pasteur pipettes

Standard 2 mL models fitting Acura® and Calibra® 2 and 5 mL macropipettes with their respective adapters. Useful alternative to macrotips, when handling PP affecting reagents.



More details

### Your local distributor

Distributed and Supported in the UK by:



24 Norman Way Industrial Estate, Over, Cambridge, CB24 5WE  
Tel: +44(0)1954 233 100 Fax: +44(0)1954 233 101  
Email: sales@camlab.co.uk Web: www.camlab.co.uk

Products and specifications are subject to change without prior notice.



 **SOCOREX**  
SWISS

Socorex Isba SA  
Chemin de Champ-Colomb 7a  
1024 Ecublens / Lausanne  
Switzerland  
Tel. +41 (0)21 651 6000  
socorex@socorex.com  
www.socorex.com

B.QTX.E – A1117