

SHEL LAB *BacBASIC*

Environmental Work Station & Incubator
For Anaerobic Applications



Ergonomic, "Gum Drop" Arm Port Doors offer maximum reach & comfort!*

***ARM PORT
DOOR DESIGN
PATENT
PENDING**



Removable viewing shield for introduction of equipment, cleaning, decontamination or maintenance.

The only thing "*basic*" about the SHEL LAB *BacBASIC* is its very reasonable price.

The *BacBASIC* is the ideal, entry level chamber for performance of truly anaerobic microbiology. Move beyond jars and the unnecessary exposure of your anaerobes to deadly oxygen.

Already have a chamber? Then the *BacBASIC* offers the perfect auxiliary chamber option for over-flow work, extra plate storage and emergency back-up for your primary chamber!

Camlab Limited
Camlab House, Norman Way
Over, Cambridge
CB24 5WE

Constant thinking.


SL
SHEL LAB

A Product of
SHELDON MFG. INC.

SHEL LAB BacBASIC

The SHEL LAB **BacBASIC** offers:

- An ample, ergonomic, user-friendly work station with a minimal, space-saving exterior footprint.
- Excellent lighting and visibility.
- A quick-entry, gloveless sleeve system.
- A rapid* Auto-Purge Pass Box for introducing Petri plates, specimens and supplies. (*Purge time less than 50 seconds.)
- Incubated storage capacity for up to 150 Petri plates.
- Incubation temperature range from ambient temperature to 42°C.
- Easy-Change anatox and catalyst packs and the catalyst can be used for up to six months before replacement or reactivation.

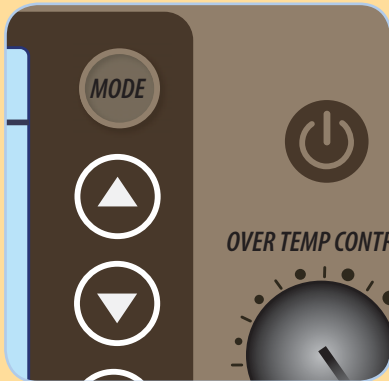

Simple Push Button Operation



Easy Change Anatox & Catalyst Packs



Central Pass Through Box



Interior Auto-Purge Button

For more information on the
SHEL LAB BacBASIC:
Call: +44 (0)1954 233100
Go To: www.camlab.co.uk

SHEL LAB **BacBASIC**

Exterior Dimensions (w x d x h)	34.75 in x 29.5 in x 29.25 in 88.3 cm x 74.9 cm x 74.3 cm
Interior Dimensions (w x d x h)	32.5 in x 19.25 in x 22 in 82.6 cm x 48.6 cm x 55.9 cm
Weight	230 lbs / 104.3 kg
Temperature Range	Ambient to 42°C
Temperature Uniformity	+/- 0.5°C